



EPILYUM
Hair Removal

Epilyum hair removal machine

Operating manual 2024



Contents

1. Safety	1
1.1 Intruduction	1
1.2 Brief introduction	2
1.2.1 Eye protection	2
1.2.2 Patient protection	3
1.2.3 Operator protection	3
1.2.4 Treatment room Requirements	3
1.3 Main warning	3
1.4 Warning labels	5
1.4.1 Warning labels on machine	5
1.4.2 Warning labels in treatment room.....	6
1.5 Laser, electrical and mechanical safety	6
1.5.1 Laser safety	6
1.5.2 Electrical and mechanical safety	7
1.6 Fire Hazards	8
1.7 Safety setting	9
1.8 Safety class	9
2 Machine introduction	11
2.1 Brief introduction	11
2.1.1 Machine outlook	11
2.1.2 Handle piece outlook	11
2.1.3 Machine specification	12
2.2 Parameter	12
2.3 Working environment	13
3 Installation.....	14
3.1 Installation environment requirements	14
3.2 Unpacking and inspection	14
3.3 Installation requirements	15
3.4 Movement	16
3.5 Transportation requirements	16
3.5.1 Original package	16
3.5.2 Other package	16
4. Operation instruction	18
4.1 Instruction	18
4.1.1 System start-up	18
4.1.2 User login	18
4.1.3 Treatment mode and parameter setting	20



4.1.4 System ready	22
4.2 Other functions.....	22
4.2.1 Track customer treatment record	22
4.2.2 Log out/Switch operation	23
4.2.3 Save or modify treatment parameter	23
4.2.4 System setting	24
4.3 Treatment	25
5. Maintenance and repair	
27 5.1 Machine cleaning	27
5.2 Water changing	27
5.3 Regular maintenance	28
5.4 Troubleshooting	29
5.5 Attention	30
6. Application guide	
31 6.1 Indications and contraindications	31
6.2 Side effects	31
6.3 Post-hair removal care	32



Chapter 1. Safety

It shows the safety considerations when using the laser hair removal system in this Chapter, especially the safety for the optical and electrical system.



Warning

It may cause personal injury if not adjust or use the machine as required by this manual. So both the operators and the technicians who use this machine must be familiar with the safety requirements and

1.1

Introduction

The EPILYUM laser hair removal system is a precise laser medical system, which is suitable for aesthetic medical applications for hair reduction treatment. There are strict safety inspections for the machines before leaving the factory, it will not cause safety problems if follow the instructions. In order to protect the operators or patients from being hurt by the laser radiation, please read this instruction carefully, and the related operators should have professional training for laser medical devices.

It is class A laser device for the EPILYUM laser machine.

The class A laser is high power laser, during usage, please pay more attention in

case of danger. Avoid any flammable items close to the laser beam absolutely.



1.2 Brief introduction

Safety is one of the main concerns of product design that should also be alert on.

Please find partial safety instructions and precautions as follows:

- ◆ The laser beam only emits from the front face of the sapphire wafer, and the energy is conducted to the skin and hair through enclosed sapphire wafer.
- ◆ The independent safety circuit will power off automatically when there is much energy.
- ◆ The Emergency Off button can cut off the power instantly
- ◆ The key cylinder switch can abstain from misuse.

1.2.1 Eye Protection

Warnings

Do not observe the laser beam or the reflected beam directly, otherwise irreversible damage may be caused to the retina.

Nominal Ocular Hazard Distance (N.O.H.D.)

This is the distance from the source at which the intensity or the energy per surface unit becomes lower than the Maximum Permissible Exposure (M.P.E.) on the cornea and on the skin. The laser beam can thus be considered as dangerous if the operator is closer from the source than the .N.O.H.D

When looking at the beam with an optical system, one has to consider the possible higher intensity entering the eye, and thus to expand the evaluated N.O.H.D.

The place to store the safety glasses should be marked with striking signs, ensuring all the personnel to put it on before powered on.

EPILYUM is equipped with safety glasses with proper wavelength.

1.2.2 Patient Protection

The safety of the patient mainly lies on a practiced beautician and a proper office Treatment. The patient should also have knowledge of the basic principles.

The patient should wear the protective glasses and know the precautions.

1.2.3 Operator Protection

The operator should read this manual carefully and receive the training of EPILYUM USER MANUAL , The operator should also know the correct operation and the possible risks.

1.2.4 Treatment room Requirements

There should be clear warning signs in the office.

When class A laser machine works, there should be warning signs on each entrance of the office.

- ◆ There should be laser safety warnings signs on each door with the wavelength

- ◆ There should be warning light on the entrance. When the machine starts, the Light is on.

- ◆ Ensure there is no highly reflective material in the office, otherwise the radiations will do harm to the body.

1.3 Main Warnings

EPILYUM can only be operated when the system inspection is done and no alarm generated. Please carefully read the following warnings

Attention!

- ◆ Laser beam radiation through window should be avoided when device is in use.
- ◆ Do not open the cover to adjust or repair the system by yourself. Only the authorized person by EPILYUM could do so.
- ◆ To make sure that rated voltage , agree with the commercial voltage (AC200-240V) in your country.
- ◆ Power off the machine and disconnect the main lead before maintaining the machine, otherwise personal injury or damage to the machine may be generated.
- ◆ Keep the treatment head clean
- ◆ Do not use the machine if the treatment head leaks; if the machine is on, please powered it off immediately.

Warnings about laser

Avoid laser beams' exposing on eyes or combustible materials, as it will cause eye damage or burning.

- ◆ It is strictly forbidden to use the laser one the same area for long time to avoid skin thermal damage.
- ◆ Operators and patients should wear laser goggles during laser treatment.
- ◆ Avoid looking at direct laser beams even wearing laser goggles.
- ◆ Hands should be away from the handle piece when turn on the machine.



1. Warning labels


1.4.1 Sample for warning labels

No. Labels

		Remark
		Name plate Information and technical information for machine and manufacturer
2		Laser warning label Remind there is laser radiation near the device.
3		Attention: Check the document along with the machine
4		Emergency Shut-down Label Press the button to shut down the machine in the case of an emergency
5		Laser Radiation and Type of Laser Beam label Under the control panel Invisible radiation

6		Located at the front of the treatment head There is laser beam here.
7	ON OFF	Key switch label Located around the key switch
8		Air switch label Located around the air switch
9		Warning of electric shock label On the back of the machine

1.4.2 Warning labels in treatment room

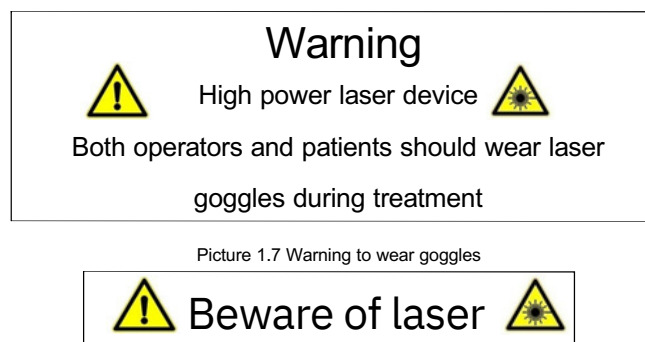

Picture 1.8 Warning for laser beam

1.5 Laser, Electrical and Mechanical Safety

1.5 Laser, electrical and mechanical safety

1.5.1 Laser safety

The Artemis EPILYUM hair removal machine is a high-tech laser device, operators should pay attention to the laser safety when use the machine, follow the below instructions:

◆Both operators and patients should wear protective goggles during laser treatment in case of eye damage. If anyone can not wear goggles, it is necessary to protect her/his eyes with some opaque materials.

◆Make sure the laser beam is closed or hard to approach.

◆Please do not turn on the machine if the case of the system is opened, or it will

cause malfunction, and make the operators exposure to laser radiation, which causes safety problems.

◆The laser system will stop emitting light to ensure the operators' safety if the environment temperature is too high, the treatment head is overheat by lack of water

or other malfunctions.

◆Only the necessary, well-trained personnel is allowed to approach the machine.

◆Make sure all staff in the treatment room know how to use and shut down the laser hair removal machine.

◆Please do not use the laser hair removal machine in the environment with explosive anesthetics or other flammable materials.

◆Please do not point the treatment head towards areas other than the target treatment area.

◆Prevent reflectors such as jewelry, watches, beauty devices or mirrors from reflecting the strong lasers.

1.5.2 Electrical and Mechanical Safety

◆Disconnect the plug after treatment.

◆Disconnect the plug before cleanse or maintenance.

◆Keep all the panels and cover plates enclosed.



- ◆Do not open the power supply or the cover plates while maintaining, or leave it untended.
- ◆The machine has is well-balanced, and can be moved with caution.
- ◆The ground electrode in the power protection sheath should be connected with earth through wall Line; A well earthed circuit is very crucial to safe operation .
- ◆Take care and to avoid scratch, while opening the outer packaging of the machine.
- ◆Check it out if the handle is well connected with the machine main body.
- ◆Check the water level before starting the machine, if it is low, add to the drain level.
- ◆Pull up the air switch and screw the key clockwise for recovery, if the machine trips because of overload.
- ◆The water should be distilled water.
- ◆Drain off the water if long time not use the machine.
- ◆Stop emitting by pressing the emergency button if necessary.
- ◆Add antifreeze to the water tank under low temperature environment.

1.6 Fire Hazards

Keep away from the inflammable materials, because the temperature will get higher once observed the laser energy.

Please clean the skin with soap and water instead of alcohol or acetone. The treatment should begin only when the skin is dry after cleansing.

Keep the laser treatment room away of papers or plastic materials, those materials will burn by absorbing enough power in certain distance.

If do not use the diode laser machine, or during the interval not use the machine,

please switch to “Standby” mode, and do not use the handle piece.

1.7 Safety setting

EPILYUM is equipped with a few safety device that the operator should be familiar with.

Key Switch: It is used to turn on or off the laser hair Removal machine. Only the key of machine can open the machine, which can prohibit illegal use. Make sure the key is pulled out after every use.

Emergency shut-down Switch: Press the red mushroom button, the machine will be shut down immediately. Slew the Emergency shut-down Switch clockwise, or the machine will keep the state of power interruption.

!

Attention: This button is only used under emergency condition to ensure safety.

Interlock: The machine can work when connecting with interlock.

Semi-automatic air switch:

1.8 Safety class

◆Protection against electric shock: Class I

◆Protection against electric shock: Type B. 

◆Laser radiation class: Class 4

◆Protection against corrosive liquid: Common

◆It cannot be used in the environment of inflammable gas.

◆Working Mode: Continuous

◆Rated voltage and frequency: Single-phase AC200-240V, 50Hz $\pm$ 10%; input power: 2000W.



The network power supply connected to this unit must use a single-phase three-wire power socket that conforms to the national standard and has good grounding .

◆There is no application part that protects against defibrillation discharge effects.

◆Non-permanent installation equipment, which is mobile.



Chapter 2 Machine introduction

This chapter is the general introduction for EPILYUM laser hair removal machine

2.1 Brief introduction

2.1.1 Machine outlook



2.1.2 Handle piece outlook



2.1.3 Machine specification

EPILYUM laser hair removal utilize the (755nm + 808nm + 1064nm) the right wavelength to penetrate into the hair follicle. Using the theory of selective laser absorption, laser can be preferentially absorbed by the melanin of the hair. And then the hair shaft and the hair follicle are heated rapidly in order to destroy the oxygen-supply organization of and around the hair follicle. The treatment window of the handle is designed with sapphire dynamic cooling system, which could bring about pain-free effect to the epidermis with pain free treatment.

2.2 Parameter

Table 2-1 Main parameter

Model	EPILYUM HYBRID ALEX SYSTEM
Wavelength	(500 - 1100)nm According to choose
Maximum Frequency	1 - 10 Hz (10 shots per second)
Laser Energy	(1 - 65) J
Device Power	5000 W
Pulse Duration	(1 - 350) ms
Spot Size	15 mm X 50 mm
Certificates	CE and ISO 9001
Handpiece Tip	Crystal Sapphire
Voltage Shot	220 V

2.3 Working environment

◆ Working environment

- A. Environment temperature : 10°C~40°C;
- B. Relative Humidity: ≤80%;
- C. Atmospheric Pressure: 860hPa~1060hPa;
- D. Power Supply: AC200~240V/50Hz.

◆ Transportation and Storage Environment

- A. Environment Temperature: 0°C~60°C; B.
- Relative Humidity: ≤100%;
- C. Atmospheric Pressure: 500hPa~1060hPa



Chapter 3. Installation

3.1 Installation environment requirements

Select a place to install the Machine

While selecting the place, you have to consider the following factors:

- ◆Keep the machine away from the inflammables, the explosives and corrosives, such as alcohol, sulfur and so on.
- ◆Keep the machine away from the heat source and the sun. Rapid change of temperature can water vapor may accumulate inside of the machine.
- ◆Put the machine in somewhere level without other goods around, and it should be at least 0.2m away from the wall.
- ◆Keep the office clean and other things in good order.
- ◆Keep the machine away from strong electromagnetic wave interference sources.
- ◆Keep away from obvious vibrations and shocks;
- ◆The office should be with good ventilation.
- ◆Keep the machine away from other goods for at least 60cm.

3.2 Unpacking and inspection

If the package is broken, please check if there is any deformation, scratch or damage of the panel. Should there be any damage, please contact us or our dealer immediately and keep the package for us to judge the reason.

Please check if the goods are included according to the following list

- ◆Main body
- ◆Handle piece



3.3 Installation Requirements

Before open the package of the laser machine, please check the whether the installation position qualified or not.

Space requirement

Electrical Requirements

The rated voltage of laser is in accordance with the requirement of the user. And the electrical requirements are as follows

Single-phase AC200~230V, 50Hz□10%

The grounded circuit contributes a lot to safe operation.

The machine should have its own transmission line and air switch. For the input power supply allows for no transient voltage, current spikes or surges.

The main safety is a 20A circuit breaker (air switch). When it is down, the power is cut off. Pull it up, the power is connected again.

Electrical interference: If use the laser machine in hospital or beauty clinics, it may cause interference to some surrounding equipment, the laser can't be used around strong magnetic fields or strong microwave fields, and should not be shared with interference devices such as MRI and CT.

Environment requirements

Air quality: There should be no corrosive particles in the air. Corrosive particles such as acids can damage the surfaces of wires, electrical components and optical components.



The dust in the air should be as less as possible. The dust can absorb light energy and be heated. The dust on the surface of sapphire may damage the sapphire when it was being hearted. The metal particles may damage electrical components.

Temperature: To ensure the best working status of the machine, the indoor temperature should be kept between 20℃ and 25℃, and relative humidity should be not exceed 80%.

Temperature: To ensure the best working condition of the machine, the indoor temperature should be kept between 20℃ and 25℃, and the relative humidity should not exceed 80%.

3.4 Movement

Please observe the following step when you need to move the machine

1. Disconnect the power cord
2. Put the treatment handle on the handle holder.
3. Loosen the brake wheels
4. Push the hand grip of the machine slowly, and do not push the display



Chapter 4. Operation Instruction

4.1 Instruction

4.1.1 System start-up

Connect to power supply, and push the air switch upwards to the position “I”, turn the key switch clockwise to turn on the machine, and go to the main interface, as picture 4.1.




Picture 4.1 Main interface






choose the air blow , you can choose one handle until the flash lamp will finish

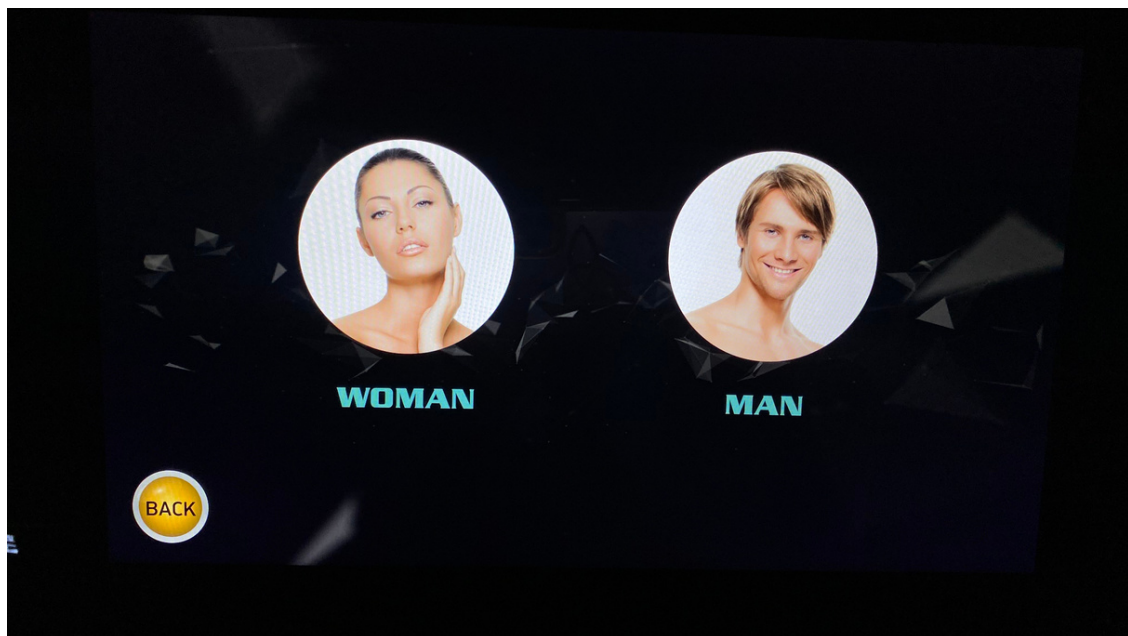

Press gender, to choose "Male" or "Female",



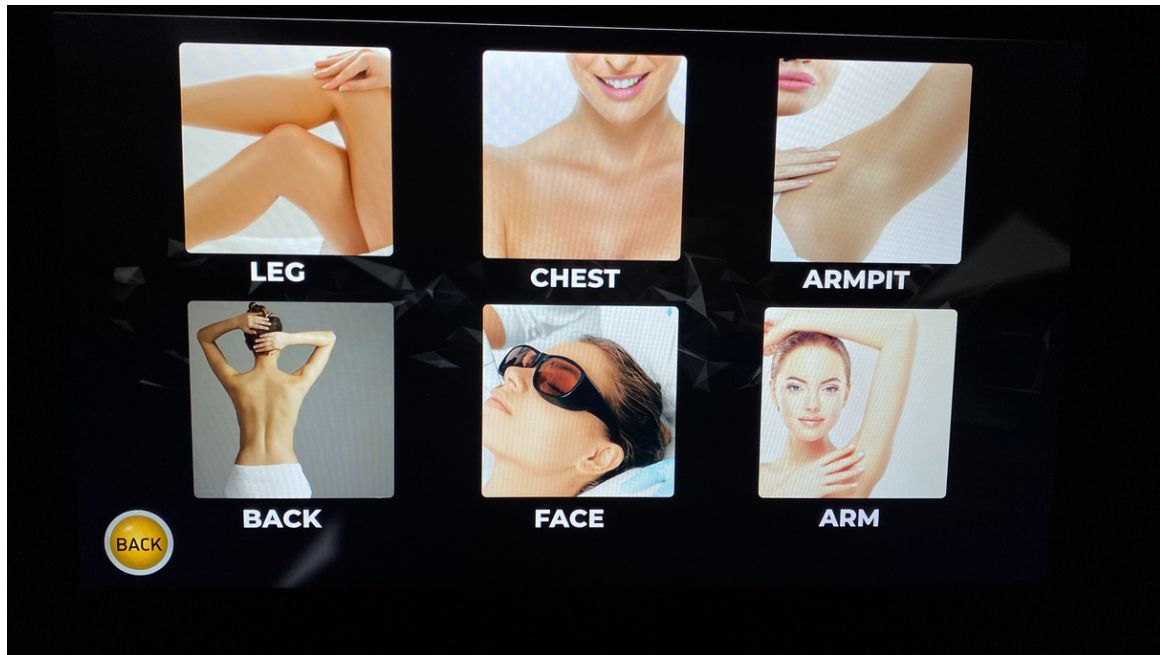


4.1.3 Treatment mode and parameter setting

Choose the treatment area, treatment mode and parameter after login the system.

- there is 6 treatment

areas for option (such as face, underarm, arm, leg, finger and bikini), choose any area

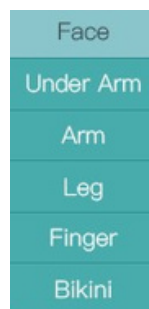


●there are 4 skin types and 3 hair types for choice. Customers can choose the right treatment mode

skin type and hair color by different treatment areas or their own requirements.

“Skin type” to choose the skin type, from darker to lighter skin.

the power density is up to 60J/cm².

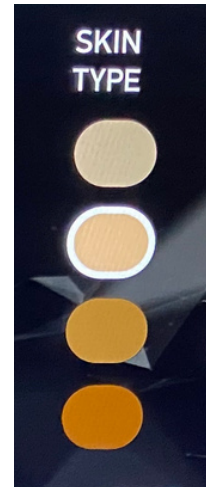

Picture 4.7 Choose skin type

Press the button “Hair type” to choose the related hair type

Picture 4.6 Choose hair type





1-frequency : its laser shots speed

choose 1 for face treatment

choose 2-3 for body treatment



2 / ENERGY - POWER / J

Tips: During the treatment, the operator should adjust the power from lower to higher, make sure the customer can tolerate, according to customer's feeling

start power in (2-3)j - in beginning flash lamp round

increase 1j in end flash lamp round


decrease 1j in darker skin area , or if you noticed pigment lesion on the skin

3 / pulse duration / ms

set pulse duration in 5ms

in case of fine hair , increase it for (5-6)ms

for air level , can press the cooling button
"" on the lower right corner.

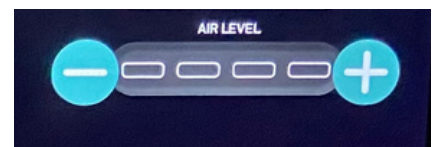


4 levels



4.1.4 System Ready

The machine is in ready status after chosen the treatment mode and set the treatment parameter, press the button “Ready” and the machine goes to “Standby” status. The handle can emit laser beam only when the machine is in Standby status. At the same time, the other setting is locked in case of disoperation.

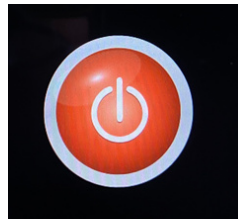

Picture 4.8 System ready Picture 4.9 System standby

Press the button “back, can go back to the search page



4.2.4 System setting


Press the button “System” you will see the system setting, it shows the machine

SN number, software version, main engine SN number, Handle piece SN number,
working time and etc. Shown as picture 4.13

4.3 Treatment

Step1: Cleanse

Clean the treatment area before laser hair removal treatment, till the target area is dry enough to do next step.

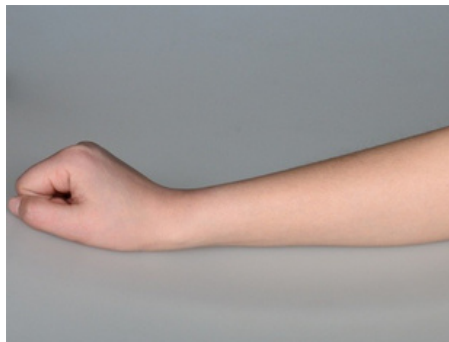

Picture 4.12 Cleanse



Step 2

:

Shave the hair before treatment.



Step 3: Set the treatment parameter

step 4 ;djust the power density and frequency.

Step 5: Hair removal

Hair removal method: The handle piece contact the skin



Chapter 5. Maintenance and repair

In this chapter the methods to maintain the machine and the solutions for common fault is included.

The user could only carry out the following maintenance work.

5.1 Machine cleaning

Clean the surface of the machine with wet cotton ball once a week. Neutral cleanser can be used, but do let the liquid soak into the inside. Clean the treatment window after each treatment.

Good habits can help you to apply the machine favorably and length the life cycle of the machine. The following suggestions are for effective clinical treatment.

Attention:

Please make sure the machine is off or in standby mode when cleaning the system.

Do not switch the machine frequently

5.2 Water changing

◆ Change the water in the water tank every 6 months to ensure the water quality,

it is suggested to use pure water.

: Procedure to change water



1. Turn off the machine and disconnect the power cord.

Check whether the machine is off and power cord is disconnected before changing water to avoid danger.

2. Lock the universal wheel

There are four universal wheels at the bottom of the machine, two of which has locks. Lock the universal wheels by treading the pedal.

3. Drain off water

Water outlet (WATER OUT) is on the back of the case, connect the water pipe to the water outlet position, and put the funnel downwards to a container, till there is no water flowing out.

4. Fill water

Connect the pipe to the water in position, and put the funnel upwards, and fill pure water to the funnel.

5. Finish changing water

Stop filling water till there is water flowing from the over flow position, and disconnect the pipe.

5.3 Regular Maintenance

To ensure the machine on its best state, periodical maintenance from professional personnel is needed, as well as daily care. It is suggest that the maintenance should be

carried out every 6 months, which including the following items.

◆Overall system inspection

◆Internal inspection, includes cleaning dust of all the components

◆Output power inspection

◆Cooling system inspection and cleaning.

◆Radiator cleaning.

◆Air fan cleaning.

◆ Inspection or replacement for the “O” shape ring of the handle connection position.

5.4 Troubleshooting

Problem	Probable causes	Solution
No display	Power supply looseness	Check whether the power cord is well connected.
	Air switch shut off	Push the air switch upward, turn on the machine again, if the air switch shut off again, please contact our after sales department.
No reflection when you press the button on the screen.	Emergency button is pressed	Contact our after sales department.
The system cannot go to main interface.	Key switch was not turned on	Contact our after sales department.
-STOP on the screen	The interlock is not connected	Contact our after sales department.
E-KEY on the screen	The temperature for the laser module is overheat	Contact our after sales department.
E-LOCK on the screen	Skin temperature overheat	Contact our after sales department.
E-LD on the screen		Contact our after sales department.
E-SKIN on the screen	Water flow alarm (Water cycle blocked)	Contact our after sales department.
E-FLOW on the screen	Water level alarm (not	Contact our after sales department.
E-STAGE on the screen		Contact our after sales department.

	enough water inside water tank)	
E-TYPE on the screen	SN number for machine main body and handle not match.	Contact our after sales department.

5.5 Attention

Attention:

1. There is high voltage when using the EPILYUM laser hair removal system, the operator should follow the user manual strictly, it is forbidden to touch the parts with electricity to avoid electric shock.
2. The system is grounded through the grounding conductor in the power cable. This protective grounding is essential for safe operation.
3. The operator should follow the user manual strictly in case of damages.
4. Do not use flammable anesthetic or oxidizing gases such as nitrous oxide and oxygen.
5. It is recommended to check the integrity of the transmission system before use



Chapter 6 Application guide

6.1 Indications and contraindications

: Scope of applications

Underarm hair removal, leg hair removal, arm hair removal and etc.

: Contraindications

Customers has below conditions can't have laser treatment such as heart disease, psychosis, hypertension, diabetes or someone who is not suitable for treatment after inspection.

6.2 Side effects

There may be side effects during treatment, and the responding is as follows;

Pain

Patients may feel different levels of discomfort, which is bearable for the adult.

The sensation of burning may last for 1h at most.

Excessive Redness and Swelling

The skin may get red and swollen after the treatment for a short period, which will be reduced in 3-7 days.

Fragility of the skin

The skin may get fragile after the treatment. Please do not make up or rub the target area.



6.3 Post-hair removal care

Cold compress rather than ice compress can release the swelling and discomfort.

The temperature of the pack for cold compress cannot be lower than 4 °C.

The blisters and ulcers on the skin after hair removal treatment can be treated with antibiotic ointments and burn ointments (e.g. Silverol or Biafine and etc). Lower concentrations of steroid ointments can also be used to reduce swelling.

Pay attention to avoid trauma within 4-5 days after hair removal treatment. Do not have hot showers, do not do aerobics or massages and etc.

If the blisters are scarred, use moisturizing ointment such as petrolatum to keep it wet.

If side effects occurred during treatment, please do not continue treatments before the target treatment area healed.

If the hair was removed partially after the treatment area healed, it is suggested to lower the power density for continuous treatments.

If there is side effects and the hair was not removed after 4 weeks, the customer should wait another 2 to 4 weeks for review. If the hair is partially removed, continuous treatment can be performed, and it is better to lower the power density during treatment.

If the hair was not removed by second review, maybe the side effect is caused by the sunburn of the patient's skin. If so, the patient should continue treatment after skin recovery.

End treatment

It is up to the beautician to decide whether to end the treatment or not. Normally when the hair was removed clearly or the patient is satisfied with the treatment result, it is OK to end the treatment.





EPILYUM PRO SERIES



CONTACT US

Next to PAR hospital , 60 m street
Business building , Erbil / Iraq

Yeşiloba Mah. 46075 Sk. E Blok
Daire No: 1 Seyhan / Adana / Turkey



+964 750 539 8000

+90 530 774 09 95



info@eslazer.com.tr

www.eslazer.com.tr



epilyumiraq@gmail.com