



PICTRON11


Picosecond Laser



- Highest Power of Any Multi-Wavelength Picosecond Laser
- Shorter Pulses (550ps) - Reduced Pain & Erythema
- Dual Pulse Width AND Four Wavelengths for Maximum Versatility
- Advanced Proprietary Technology & Elegant Ergonomics




EPILYUM PICTRON11 is the next-generation product following the top-of-the-line series, **Pictron11** offered four wavelengths: **532 nm, 585 nm, 650 nm, and 1064 nm**. This variety of options gave it the versatility to treat even the most difficult cases. It has inherited this novel feature, but operating in the powerful and stable picosecond domain. The availability of specialized handpieces at four different wavelengths, combined with the ps-domain pulsewidth, makes the removal of both multicolored and black tattoos even more effective. Moreover, it can address a variety of epidermal and dermal pigmented lesions and can also treat non-pigmented skin conditions and lesions such as wrinkles, acne scars, and enlarged pores.

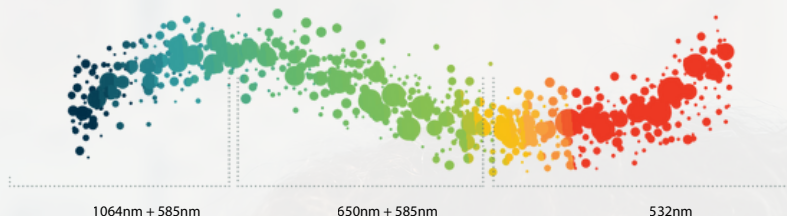

ULTRA-SHORT PULSE PERFORMANCE

PICTRON11 delivers cutting-edge precision with an exceptionally brief pulse duration of just 550-750 picoseconds—among the fastest in its class. This advanced, versatile system empowers practitioners to customize treatments by adjusting key parameters such as wavelength, energy output, spot dimension, emission mode, and repetition frequency, ensuring truly personalized results for every patient.


Pictron11 delivers a full range of laser wavelengths to cover the entire absorption spectrum of melanin, oxyhemoglobin, and water.

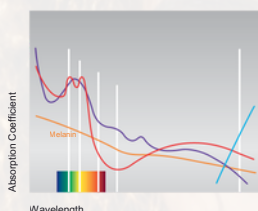

DUAL-SEQUENCE EMISSION

PICTRON11 incorporates an innovative dual-pulse delivery system that emits two precisely timed pulses in each laser cycle. This unique approach amplifies traditional picosecond laser advantages by intensifying the photoacoustic effect while introducing a controlled thermal component, resulting in improved fragmentation and enhanced overall treatment efficacy.

PICTRON11 works by sending extremely short picosecond pulses into the skin. These very fast pulses create strong energy bursts that break tattoo pigments into tiny pieces. The shorter the pulse, the finer the pigment is broken—from larger bits down to fine dust. For colored tattoos, the special double-pulse mode helps break the ink even better, using more power and stronger shock effects than with a single pulse.

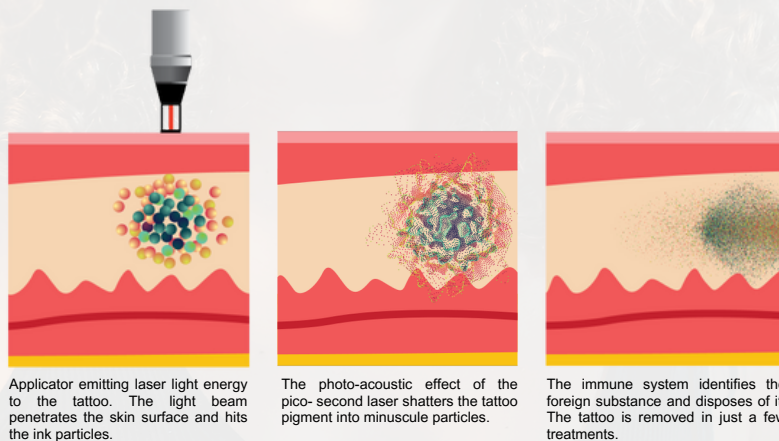

Pigmented Lesions

PICTRON11 can treat both shallow and deep pigmented spots. With its multi wavelengths. Its ultra-short pulse duration delivers very high peak energy that heats the pigment quickly, causing it to break apart through a strong photoacoustic effect. This process breaks the pigment into very small pieces, often needing fewer treatment sessions than other lasers.

WIDE VARIETY OF APPLICATORS

Pictron11 features a comprehensive range of multiple handpieces, each specifically designed for targeted treatments and clinical indications—together delivering a complete spectrum of advanced solutions.



Focus Handle (Standard handle)

Melasma, freckles, sunspots, age spots, Ota's nevus, black face, blue and black tattoos, black eyebrow tattoos, eyeliner, skin rejuvenation

Honeycomb Handle (Optional)

Skin rejuvenation, dermal spots, Ota's nevus, birthmarks

Gaussian Handle (Optional)

Freckles, sunspots, age spots, Ota's nevus, black face, blue, black tattoos, black eyebrow tattoos, eyeliner, skin rejuvenation

650nm Handle: (Optional)

650nm is used on green tattoos.

585nm Handle: (Optional)


585nm is used on red, purples, multicolor tattoos.

Fractional Handle : (Optional)


Pictron 11 features both standard and fractional handpieces, giving clinicians versatile tools for a wide variety of treatments. Its fractional mode delivers concentrated photomechanical energy in precise micro-arrays, ensuring targeted results while preserving the safety of surrounding tissue

Seborrheic keratosis, wrinkles, sagging skin, scars, acne scars, and pores.


- Enhanced generation of photoacoustic shockwaves
- Reduced frosting and pressure on surrounding tissue
- Direct pigment removal via ablation
- Better and quicker healing of fractionated skin

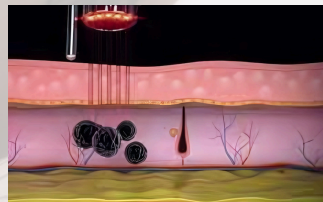

Treatment with a Fractional handpiece of PICTRON 11

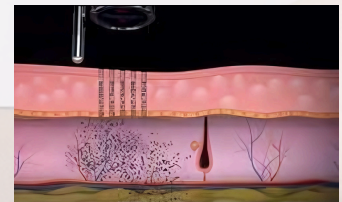

Removal of a pigment via ablation

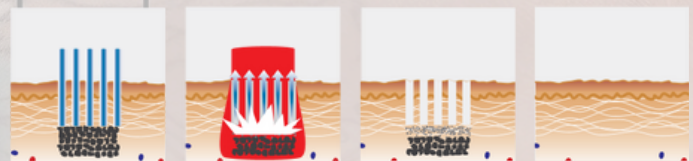

BENEFITS



Shortest Pulse Duration,
Enabling High Peak Power



Unique Double Pulse Mode



Suitable for All Skin Types



Improved ROI with Excellent
Results in Fewer Treatments



Versatile Platform Enabling Multiple
Patient-Tailored Treatment Options




offers 9 different spot sizes from
2 mm to 10 mm, making it easy to
use



FULL RANGE OF INDICATIONS

- Skin Rejuvenation
- Dark Ink Tattoos
- Colored Ink Tattoos
- Lentigines
- Acne Scars
- Wrinkles Reduction
- Pores & Texture Improvement
- Café-au-lait Birthmarks
- Benign Pigmented Lesions

Epilyum Pictron11 Treatment Progress

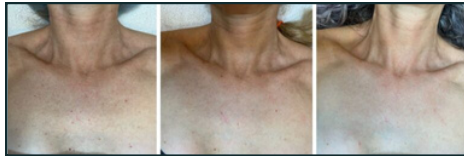

Wavelength = 532nm
Frequency = 1hz
spotsize = 2mm
power = 125mj



Wavelength = 1064nm
spotsize = 1mm
Frequency = 1hz
power = 600mj



Wavelength = 1064nm
spotsize = 2mm
Frequency = 3hz
power = 500mj



Wavelength = 532nm
Frequency = 1hz
spotsize = 3mm
power = 175mj



Wavelength = 532nm
Frequency = 1hz
spotsize = 3mm
power = 150mj



Wavelength = 1064nm
Frequency = 3hz
spotsize = 3mm
power = 600mj

Operationally Flexible and Portable

Quick readiness: Complete warm-up in less than 60 seconds

Efficient use of time and power: Turn off between sessions without needing to warm up again

Easy mobility: Operates on standard 220-volt power for use in multiple rooms

beam of light transmission: 7 articular-am of light guide transmitting

System Specifications

Wavelengths: 1064nm, 650nm, 585nm, 532nm

Spot size: 2-10mm

Pulse energy: 50-2000mj pulse

Pulse duration: 550-750ps

Pulse shape: single pulse , double pulse, continuous pulse

Frequency: Upto10 Hz

Cooling manner: Closed-off water circulation

Operation Interface: 15.6" TFT True Color LCD Screen Picosecond

Warnings and Risks: The EPILYUM PECTRON 11 picosecond laser (wavelengths depending on the applicator) is intended exclusively for use by physicians trained in the operation of picosecond laser technology. Improper treatment settings or misuse can result in serious injury to patients or operating personnel. Risks associated with any picosecond laser procedure, particularly when using the fractional mode, may include transient erythema, edema, hypo- or hyperpigmentation, infection, skin induration, blistering, or scarring. Before treatment, carefully review the EPILYUM PECTRON 11 user manual and accessory guides for a full list of intended uses, contraindications, and precautions.



For more information:
www.epilyum.com



Epilyum Ltd.
Yesiloba District, 46075. Street, E Block Flat: 1
Seyhan Adana
Mobile: [+90 539 649 20 96](tel:+905396492096)