



AxisOne CO₂



**Super Pulse
CO₂ Laser with Multi Mode
CO₂ FRACTIONAL LASER**

Fractional Mode / Ultra Mode / CW Mode

Deep Resurfacing

Rejuvenation Laser

Gynecological Applications



Cutting-edge advancements are on the horizon, ready to reshape our entire workflow.

CO₂ fractional laser is a far-infrared laser that delivers high-energy pulses focused on microscopic treatment zones. The intense heat instantly vaporizes targeted skin columns, creating precise microthermal ablation channels while inducing controlled thermal diffusion to surrounding tissue. This thermal effect coagulates collagen around the ablation zone, forming a protective barrier that reduces infection risk. Additionally, heat diffusion stimulates fibroblast activity, promoting neocollagenesis and accelerating skin remodeling.

Why Professionals Trust **AxisOne CO₂** ?



The **Epilyum AxisOne CO₂** uses advanced **RF Metal Tube** technology to provide high power output, consistent energy delivery, and superior safety—ensuring precise, reliable results for every treatment.


AxisOne CO₂ employs advanced, physician-focused technology to deliver deeper penetration with ultra-short pulse durations, significantly reducing thermal damage to surrounding tissue. This ensures highly effective treatments with improved safety profiles and superior aesthetic results.

Epilyum is an advanced, physician-tested brand committed to delivering exceptional quality and clinical excellence. Through rigorous testing and continuous innovation, Epilyum ensures that every device meets the highest medical standards for safety, reliability, and effectiveness. Trusted by doctors worldwide, Epilyum empowers practitioners to deliver superior results with complete confidence.

Super Pulse CO₂ gives you the ability to treat both superficial and deep skin layers simultaneously with precision control over the intensity, pattern, and depth of ablation. Both function can work at same time and independently. All was controlled by operator, which is easy for treatment.


Fast speed

Effectively improve the speed of random filling of graphics and greatly reduce the pain of treatment.

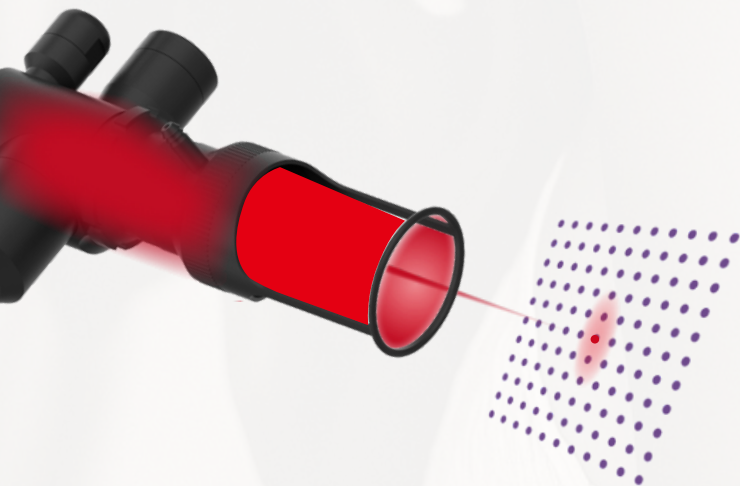


Special SuperUltra

delivers high peak power in a short pulse duration for deep penetration with ideal ablation and coagulation ratio. It make laser fully power in short time, total power will be add 20% more in same setting.

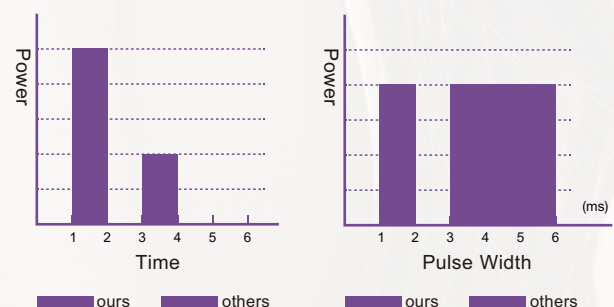

Different Graphics Modes


Deep Mode : Delivers energy deep into dermal tissue, achieving a powerful fractional effect through the precise combination of controlled ablation and effective coagulation for enhanced remodeling.

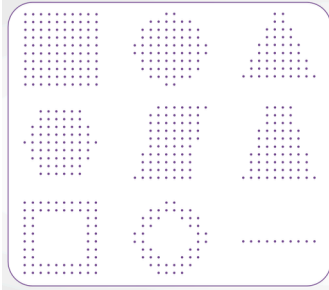


Ablative Mode : Delivers precise superficial fractional treatments in the epidermis to stimulate cellular renewal, while also penetrating deeper to achieve an intensified fractional ablative effect on structures below the dermal-epidermal junction, making it ideal for delicate areas such as eyelid tightening.




Combo Mode: Utilizes a sophisticated scan pattern that synchronizes varied penetration depths with controlled thermal energy, delivering optimized, simultaneous treatment of both epidermal and dermal layers for superior clinical results.

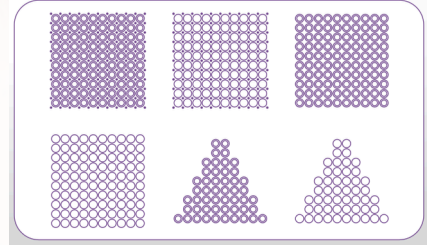


Multiple Filling Modes

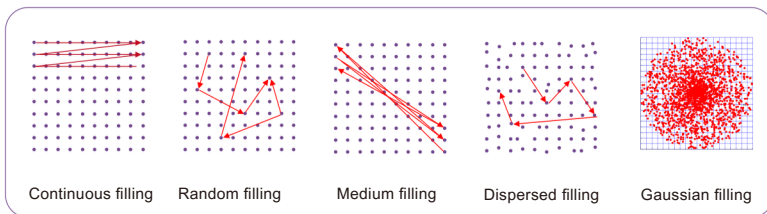

Random Mode : Employs randomized fractional channels to provide powerful, low-pain therapy, maximizing efficacy while enhancing patient comfort.

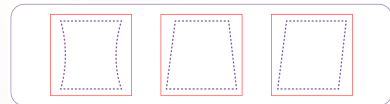

6 KINDS OF TREATMENT MODE

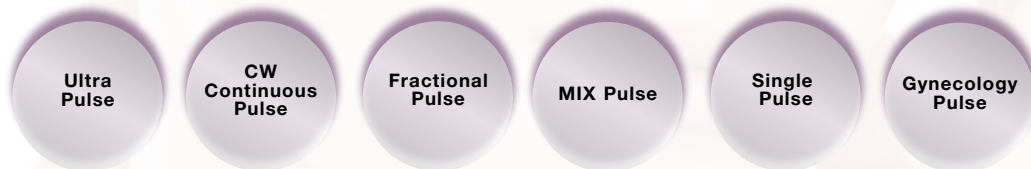

■ Fractional Handpiece → Fractional Resurfacing Module



Ultra Pulse : When ablation depth is equal, more energy means more excess heat in the tissue. **AxisOne CO₂** UltraPulse technology enables the deepest penetration with the lowest energy. This reduces patient discomfort, downtime and thermal damage. **AxisOne CO₂** also enables the shortest pulses per energy level further enhancing patient comfort.

Skin tightening, Wrinkle, Large pore, Acne scar, Skin resurfacing, Stretch mark, Laser peeling, Rubbing mode, Pudenda, LVR.

■ Surgical Handpiece → Precision Surgical Module



The **AxisOne CO₂** laser is also equipped with a dedicated surgical handpiece designed for precise, controlled tissue incision, enabling advanced cutting applications with minimal thermal damage.

Corns, Warts, Freckles, Birthmarks(mole), Phimosiectomy, Hemorrhoids and anal fistulas, Condyloma, Snoring, Skin resurfacing, Tumor elimination, Fine wrinkle, Syringoma.

■ Gynecology Handpiece → Gynecological Treatment Module



It supports a broad range of gynecological treatments, enabling clinicians to manage diverse conditions with advanced accuracy and patient safety.

Excision of CIN and VIN, vaporization of genital warts, fractional treatment for GSM and laxity, remodeling lichen sclerosus, scar revision, cervical stenosis treatment, intraoperative hemostasis.

■ protective tip - reducing the risk of deep-tissue injury.



Designed for subcutaneous placement, this protective tip mitigates thermal transmission to deeper tissues, enhancing treatment safety.

Treatment Outcomes :



WAVELENGTH

1600NM

LASER TYPE

CO₂ LASER RF TUBE (Ultra Pulse, CW, Fractional Pulse, MIX Pulse, Single Pulse, Gynecology Pulse)

DEVICE POWER

1-60 W

SCAN RANGE

MAX 20mm * 20mm

SPOT DIAMETER

100um - 2000um

SCAN MODE

ORDER, CENTER , RANDOM

TREATMENT MODE

Fractional mode, Pulse mode, Gynecology mode

OUTPUT GRAPHICS

Circle, ellipse, triangle, rectangle, square, hexagon, straight line

BEAM DELIVERY

7-joint articulated arm with rotation up to 360°

VOLTAGE

220/110 V, 50/60 HZ

SMOKE EXHAUSTING SYSTEM

Air Cooling

Warnings and risks : CO₂ lasers are intended solely for use by physicians trained in the use of the Carbon Dioxide laser (10.6 μm) wavelength. Incorrect treatment settings or misuse of the technology can present risk of serious injury to patient and operating personnel. Risks that may be associated with any CO₂ laser procedure may include change of pigmentation, infection, erythema, skin induration or scarring. Read and understand the CO₂ systems and accessories operator manuals for a complete list of intended use, contraindications and risks. The use of EPILYUM AxisOne CO₂ laser is contraindicated where a patient has taken Accutane (Isotretinoin) within the past 6-12 months, has a history of keloid formation and demonstrate excessive or unusually prolonged erythema

