

PRO SERIES 2025

**EXPERIENCE THE CONFIDENCE OF
FLAWLESS SKIN WITH OUR
CLINICALLY PROVEN LASER
TREATMENTS**




Brand Values

Epilyum is committed to science, proven results, and trust, ensuring clinically validated outcomes that strengthen our long-term partnerships with physicians, their clinics, and the patients they serve.

It Combines


Epilyum sets a new standard in laser hair removal with unmatched speed, power, and precision, delivering superior efficacy, seamless usability, and the highest levels of safety and patient satisfaction.

- ✓ **Advanced Air-Cooling System (-35 C)** ❄️
- ✓ **Dual-Hand Technology**
- ✓ **Frequency (3hz) , Fullbody time = 30 min** ⚡️
- ✓ **Rectangler Spot Size Design**
- ✓ **Inverter Motor , Energy Saving (22)A** ⚡️
- ✓ **Silent System**
- ✓ **Two Silicone spot size cover (Face + Body)**
- ✓ **Broadband wavelength (500nm - 1200nm)**
- ✓ **Low Cost & High Profit , Return on Investment (ROI)**

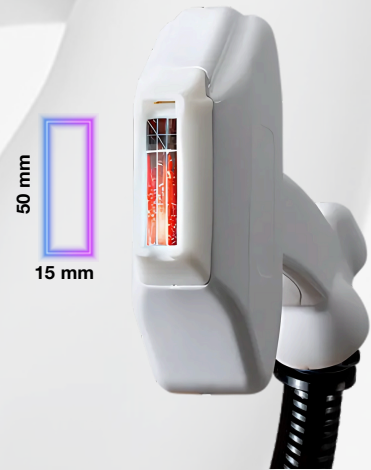

Revolutionary Rectangle Spot Size

The revolutionary rectangle spot provides faster and uniform coverage of the treatment area, with larger spots , penetrating deeper. no overlap, no burns and no retouching missed areas .

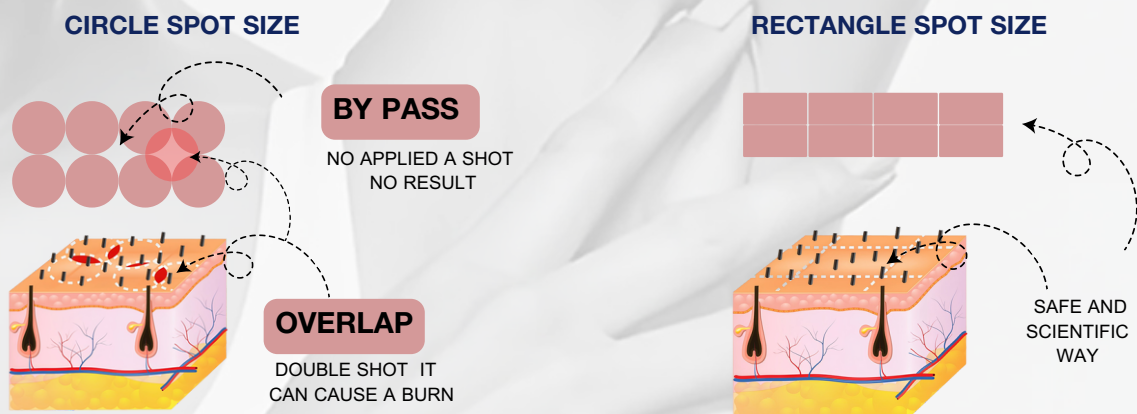

★ **9**
YEARS OF EXPERTISE

👤 **8.000 +**
10K PLUS CUSTOMER
SATISFACTION

🤝 **%99**
POSITIVE FEEDBACK

DUAL-HAND

Epilyum features a dual-handpiece design that minimizes service downtime, ensuring continuous operation and maximum productivity.



Advanced AIR-COOLING System

Epilyum integrates advanced air cooling technology to deliver continuous epidermal protection and superior patient comfort, enabling safe, high-efficiency treatments.



SPOT COVER

Hygiene process

The Epilyum laser features anti-bacteria silicone spot cover, a single-use component ensuring hygiene and optimal treatment distance between the lens and the skin.

Full-Body Spot Cover :

Designed for larger areas, maximizing efficiency.



Facial Spot Cover :

Smaller for precise targeting in delicate areas.



Advanced New IPL System for Comprehensive Treatment

The upgraded Epilyum IPL system (500–1200nm) is designed for highly effective hair removal while addressing various skin concerns with precision. It ensures optimal results across all skin tones and hair types by targeting hair follicles at every growth phase for long-lasting effectiveness.



80% High Power & Broad Spectrum Coverage

With fluency up to 100 Joules, Epilyum effectively treats various skin types and hair colors, achieving up to 80% hair reduction after the first session, while subsequent sessions serve as maintenance treatments


Time-Efficient

Epilyum's advanced technology, with a large spot size and high frequency rate, completes full-body laser hair removal in just 30-40 minutes—fast, effective, and comfortable for smooth, hair-free skin.





Cost-Effective

The Epilyum laser is highly economical for doctors due to its reduced operational costs. The primary cost is the flash lamps, making it a highly profitable investment. Doctors can enjoy up to 95% profit margin, with only 5% of the income allocated to the device's cost. This makes Epilyum a cost-effective choice for practitioners.

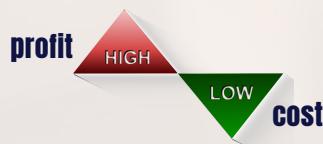


High-Definition Microscopic image synthesis technology

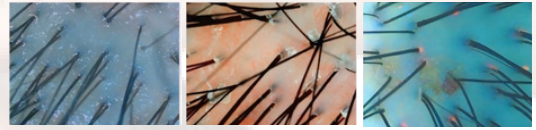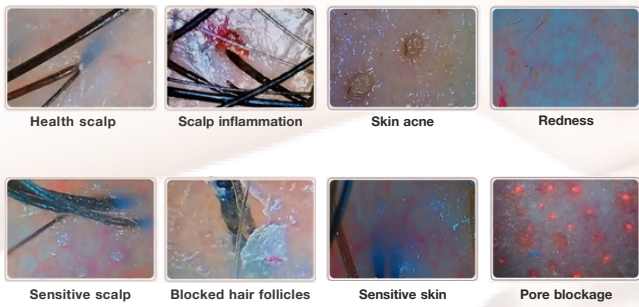

- The skin moisture analyzer can be used for detecting skin problems such as pimples, color spots.
- Suitable for testing scalp oil content, hair moisture, flexibility, etc.
- Good for solving scalp and hair follicle problems.


Superior Effectiveness on Fine Hair

The effectiveness of fine hair removal is largely influenced by pulse duration, as illustrated in the provided diagrams.

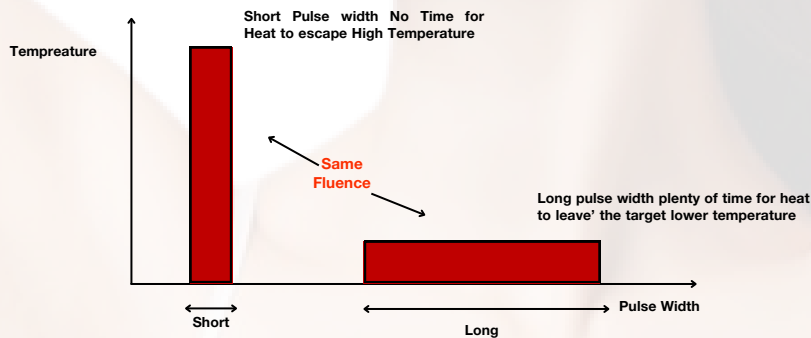


1. Shorter Pulse Durations

higher peak power and increased temperature. This rapid heat buildup prevents heat from escaping, ensuring effective follicle destruction in fewer sessions.

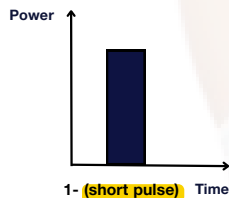


2. Longer Pulse Durations

more heat dissipation, reducing the peak temperature on the skin surface. While this is safer for darker skin types, the energy is spread out over a longer period, leading to insufficient thermal damage for effective fine hair removal.

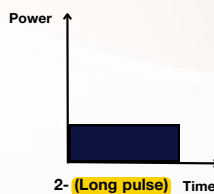

3. Pulse Train Duration

Optimization in IPL Systems Pulse trains (multiple short pulses with delays) can enhance safety while still maintaining effectiveness.

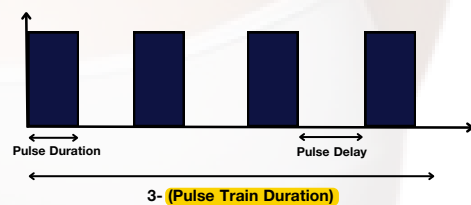


Direct Result after laser treatment



Pulse duration = 6ms
Frequency = 1hz
power = 4J



Pulse duration = 6ms
Frequency = 1hz
power = 4J

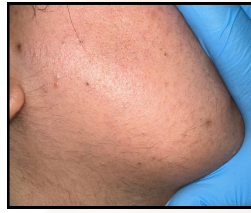

Pulse duration = 6ms
Frequency = 1hz
power = 4J

TECHNICAL FEATURES

MODEL

EPILYUM PRO SERIES 2025

WAVELENGTH

(500 - 1200)nm
According to choose

MAXIMUM FREQUENCY

1 - 10 Hz (10 shots per second)

LASER ENERGY

(1 - OVER 100) J

DEVICE POWER

5000W

PULSE DURATION

(1 - 8) ms

SPOT SIZE

15 mm X 50 mm

CERTIFICATES

CE and BELTCERT

VOLTAGE SHOT

220 V



+964 750 539 8000
+964 774 405 2020



www.epilyum.com
info@eslazer.com.tr
www.eslazer.com.tr



info@epilyum.com



epilyum.iraq_



epilyum_iraq