



EPILYUM REVIVE

INTENSE PULSE LIGHT 2026



Epilyum IPL

Epilyum IPL is an advanced Intense Pulsed Light (IPL) system designed to address over **20 skin concerns** including acne treatment, pigmentation correction, vascular lesion removal, skin rejuvenation, and hair removal. Utilizing Optimal Pulse Technology (OPT) and Multiple Sequential Pulsing (MSP), Additionally with its advanced cooling system Epilyum ensures safe, effective, and precise treatments with minimal discomfort.

Treatment Indications & Wavelength Specifications

1. Acne & Oil Control

- Wavelength Range: **420-600 nm**
- Targets *Propionibacterium acnes* (P. acnes) bacteria, reduces sebaceous gland activity, and controls inflammation, leading to the Reduction in acne lesions, oil control, and smoother skin texture.

2. Pigmentation Correction

Epidermal Pigmentation (Sunspots, Freckles, Age Spots)

- Wavelength Range: **515-755 nm**
- Selective absorption of melanin for pigment fragmentation and natural skin renewal, leading to the lightening of dark spots and even skin tone.

Dermal Pigmentation (Melasma, Post-Inflammatory Hyper pigmentation)

- Wavelength Range: **550-755 nm**
- Deeper penetration to break down stubborn melanin deposits, leading to the Gradual reduction in deeper pigmentation concerns.

3. Vascular Lesions Treatment

Rosacea, Flushing, and Fine Vessels

- Wavelength Range: **530-590 nm**
- Absorbed by oxyhemoglobin, causing coagulation and gradual clearance of visible vessels, leading to the reduction in redness and vascular lesions.



Larger Facial and Leg Veins

- Wavelength Range: **590-1200 nm**
- Deeper penetration for larger blood vessel coagulation, Resulting in the Fading and elimination of prominent vascular lesions.



4. Skin Rejuvenation & Anti-Aging

- Wavelength Range: **560-1200 nm**
- Stimulates collagen production, tightens skin, and improves texture, leading to Reduction in fine lines, wrinkles, and pore size.




ADVANTAGE

Epilyum IPL delivers high-precision, professional-grade results for dermatologists and aesthetic practitioners. Its advanced filter technology, cooling system, and customizable settings make it an ideal choice for effective, safe, and comfortable IPL treatments.


OPT and MSP Technology

OPT - Optimal Pulse Technology

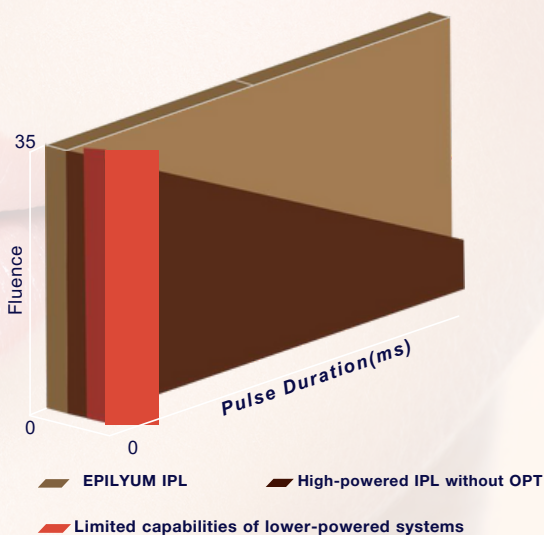

- High peak power, shorter pulses - optimal for IPL skin treatments using photo rejuvenation and the treatment of benign pigmented lesions.
- Homogeneous power over long pulse durations.

MSP - Multiple-Sequential Pulsing

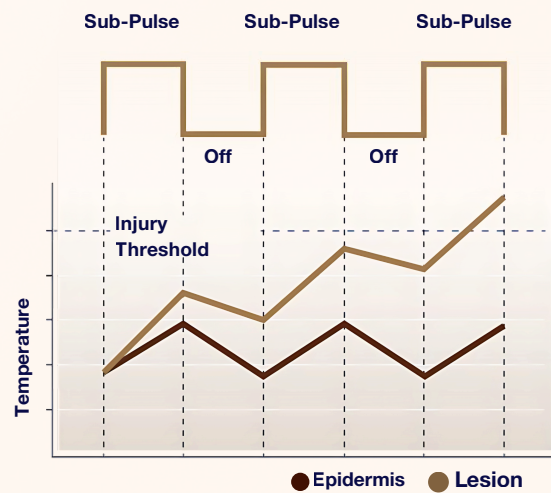


Multiple-Sequential Pulsing allows for cooling the skin in between pulses, reducing the chance of skin damage during treatments.

Quick-Change Spot Size Adaptors

Quickly and easily change spot sizes with magnetic gold pieces that instantly adjust for comfortable treatment of both small and highly contoured areas.

IPL Filters and Clinical Applications

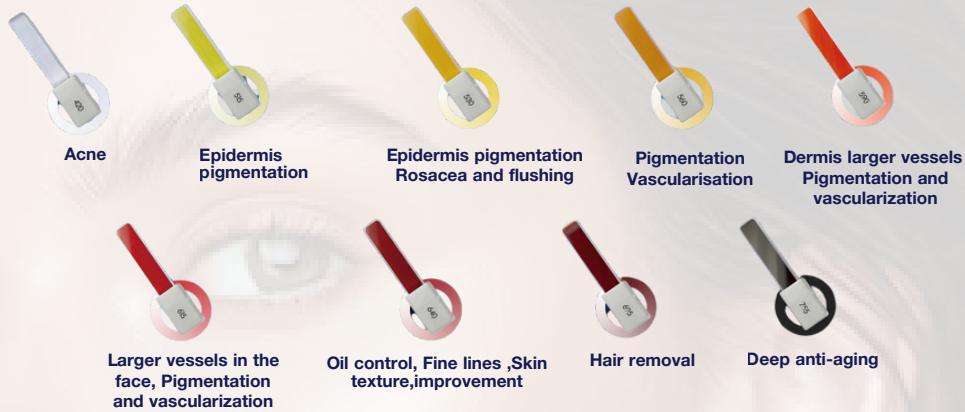

TECHNICAL FEATURES

MODEL	Epilyum Revive 2026
WAVELENGTH	(420 - 755)nm According to choose
FLUENCE	5-25 J/cm ²
SPOT SIZE	10*40 mm
INPUT POWER	3000W
PULSE DURATION	(1 - 350) ms
CERTIFICATES	CE and ISO 9001
VOLTAGE SHOT	220 V
DEVICE COOLING TEMPERATURE	-4-20°C Adjustable

Warnings and Risks: IPL (Intense Pulsed Light) systems are intended solely for use by qualified dermatologists trained in IPL technology. Incorrect treatment settings or misuse of the machine may present a risk of serious injury to both the patient and the operator. Potential risks associated with any IPL procedure may include skin burns, post-inflammatory hyperpigmentation (PIH), post-inflammatory erythema (PIE), changes in pigmentation, infection, redness (erythema), swelling (induration), or scarring. Always read and fully understand the IPL system and accessory operator manuals for the complete list of indications, contraindications, and possible risks.

EPILYUM
DISCOVER THE CONVENIENCE



For more information:
www.epilyum.com



Epilyum Ltd.
Yesiloba District, 46075. Street, E Block Flat: 1 Seyhan
Adana