



ThermiLif HIFU

Three Technologies, One Machine

Bipolar RF

SD - synergy dotting

MFU - Micro Focused Ultrasound



ThermiLif HIFU

is a machine integrates innovative, world-first technology to deliver more focused and effective treatments. By combining two powerful energy sources, it achieves exceptional results in skin lifting and shaping.

✓ MFU & RF Synergy


The world's first technology combining **Bipolar RF, SD (Synergy Dotting)**, and **MFU (Micro-Focused Ultrasound)** delivers a more powerful, precise, and comprehensive rejuvenation treatment, ensuring optimal results across multiple skin layers for lifting, tightening, and collagen regeneration.

✓ precision Lifting

Utilizes a point focusing lifting technology for highly targeted energy delivery to achieve superior lifting and tightening.

- ✓ The integrated **Bipolar Radiofrequency (RF)** system induces controlled dermal heating to stimulate neocollagenesis and elastin remodeling, resulting in progressive skin tightening and rejuvenation. Simultaneously, the **Micro-Focused Ultrasound (MFU)** technology delivers precise ultrasound energy to targeted depths, creating thermal coagulation points and micro-necrotic zones within the SMAS and muscular layers.

This dual-modality approach enhances tissue contraction and collagen synthesis at multiple depths, offering **superior lifting and tightening** outcomes compared to conventional single-technology systems.

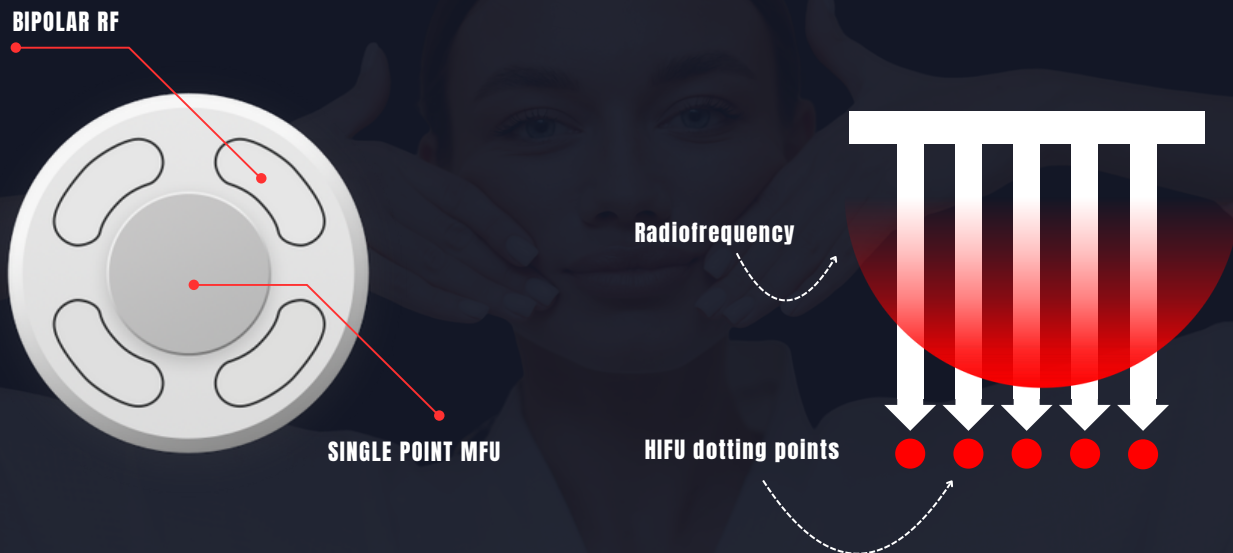


2-Layer Safe Operation Guarantee

Smart SD Rotation Sensor:

An intelligent sensor is built in to prevent repeated treatment of the same area, ensuring a safer procedure.

Smart Contact Sensor:

The machine includes a sensor that ensures energy is only output safely, likely when full contact is made.

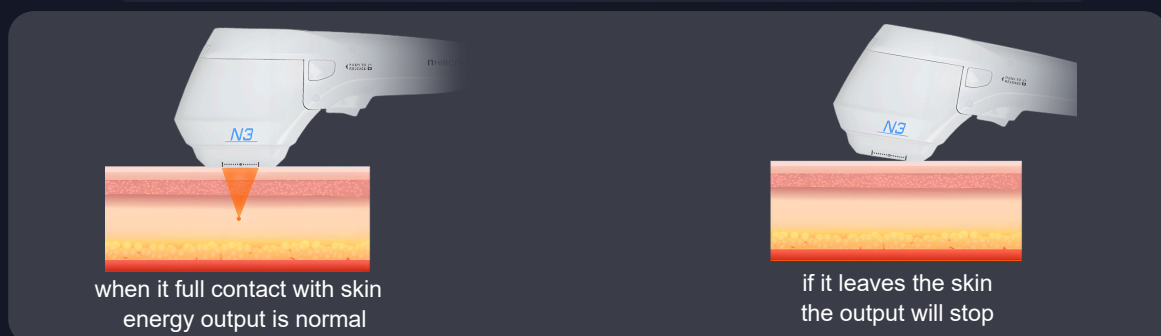

4-in-1 Versatile Applications & Precision Cartridges

ThermiLif HIFU 4-in-1 provides unmatched versatility with a range of advanced cartridges and features designed for customized, effective treatments. The system is equipped with 4 sets of therapy sticks and diversified treatment heads that utilize new technology.

- ✓ Reduces change-over time by allowing quick switching between the 4 included treatment heads.
- ✓ Features new MFU & RF synergy technology for easy, safe treatment of the face, neck, and body.

SD
4 Types



Adjustable Depth Control:


Practitioners can use the adjustable depth feature to precisely target different layers of tissue for optimal results.

FL
6 Types



Progressive Lifting and Multi-Layer Targeting

Instant result reach to 50%, after treatment result reach to **70%** (peak result will be in the third-sixth month after treatment). The depth can reach deep fascia SMAS layer can treat the epidermis, dermis, fascia and adipose.

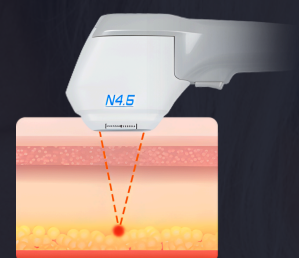




LCD Display:

Simple screen for setting and adjusting energy output, with touch switch control.



Reduced Treatment Time:

Simultaneous MFU and RF treatment with easy head switching speeds up the process.

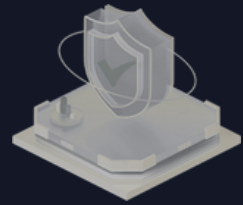

Tailor-Made Treatments:

Pre-calculated dose parameters (time/energy) ensure quick, correct, and safe settings.



CARTRIDGE KIND 9mm
ENERGY (j) 1.4j
NUMBER OF SHOTS 800 lineShot
4 SESSIONS



CARTRIDGE KIND 4.5mm
ENERGY (j) 1j
NUMBER OF SHOTS 1000 DoteShot
2 SESSIONS



CARTRIDGE KIND 4.5mm
ENERGY (j) 1.5j
NUMBER OF SHOTS 110 LineShot
4 SESSIONS



CARTRIDGE KIND 2mm
ENERGY (j) 1.2j
NUMBER OF SHOTS 100 lineShot
3 SESSIONS

TECHNICAL FEATURES

Technology	HIFU 4 IN 1 Machine
Operation Modes	Dotting mode , Fractional mode
Treatment Heads	4 sets of therapy sticks
Depth	Adjustable depth to treat the target layer
User Interface	LCD display
Safety Certification	2-layer safe operation guarantee certification
Security Feature 1	Dual patent approved security
Security Feature 2	SD rotation sensor
Primary Functions	Remove wrinkles ,Lifting and tightening skin
Line cartridge depth	1.5mm, 2.0mm, 3.0mm, 4.5mm, 6.0mm, 9.0mm
Point cartridge depth	1.5mm,3.0mm,4.5mm, 9.0mm

Warnings and Risks: HIFU (High-Intensity Focused Ultrasound) devices are intended solely for use by qualified dermatologists. Incorrect settings or misuse may cause serious injury. Possible side effects include redness, swelling, tenderness, bruising, numbness, or tingling. In rare cases, burns, fat atrophy, or nerve injury (including facial nerve damage) may occur due to improper energy delivery. Always follow the HIFU system manual and safety guidelines carefully.

