



AI REVEAL

ARTIFICIAL
INTELLIGENCE
SKIN ANALYZER



EPILYUM AI REVEAL SKIN ANALYZER

In today's fast-paced world with a research and development focus on solving skin problems, **Epilyum AI Reveal** integrates **eight** spectral imaging technologies and can professionally and objectively analyze **seventeen** problems of facial skin with flexible operations. The original intention of the research and development is to take photos and analyze reports with just one click, making it more convenient to operate.

Professional Skin Analysis:

Harnessing the power of 8-spectrum imaging technologies and intelligent AI algorithms, **EPILYUM AI REVEAL** provides a comprehensive, data-driven analysis of your skin. It identifies up to **17+ skin concerns** including wrinkles, acne, spots, pores, pigmentation, and more—instantly and accurately.

Visual Communication & Deep Diagnostics

The device generates high-resolution **3D facial imaging**, mapping volumetric changes, skin tone, and texture. This creates a shared visual language between practitioners and patients, enhancing communication and treatment trust.

Personalized treatment Guidance

- Individualized skincare plans based on AI analysis
- Precise product recommendations tailored to your skin's unique condition
- Dynamic monitoring of skincare progress over time

Recommended Treatment Products:

Based on your skin analysis, our system recommends the most suitable treatment and skincare products, ensuring optimal effectiveness and quality. Designed for ease-of-use, **EPILYUM AI REVEAL** offers one-click photo capture and analysis, multilingual support, and intuitive interface navigation—ideal for both clinics and beauty professionals.

Seventeen testing indicators



Eight spectral images






Analysis of 15 major symptoms and more than 30 detection dimensions


Match the corresponding pattern diagram. Beginners only need to analyze the sub-healthy dimension based on the trend chart.

Analysis of aging Sensitive analysis Pigment analysis Skin quality analysis

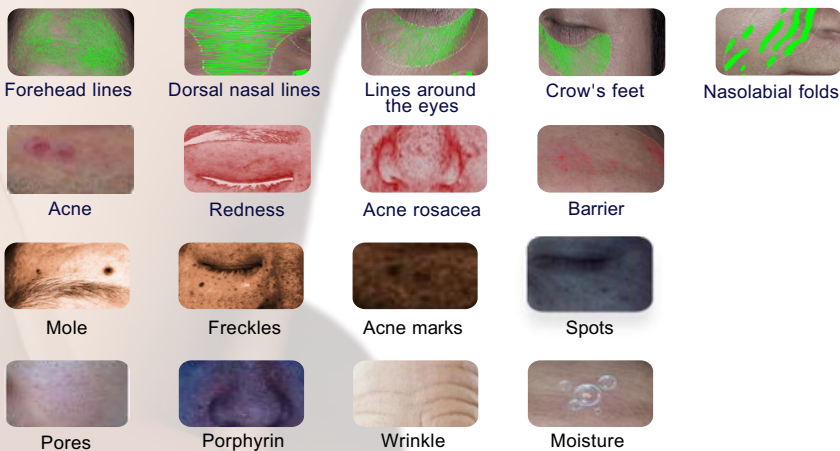


Clinical Skin Concerns :

Sebum

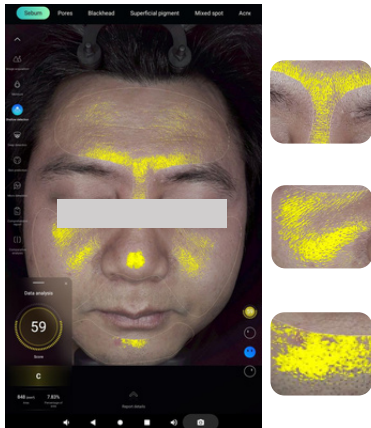


The oil secretion of the skin surface can be checked under positive polarized light source. The algorithm displays areas of the skin with active oil secretion through yellow fluorescence, through the form of data, you can see the oiliness of facial skin more clearly and intuitively.

PORES

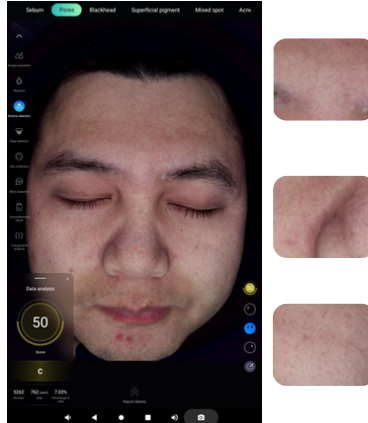


Under negative polarized light source, it is possible to check if enlarged pores have formed on the skin surface. The algorithm uses RBX technology to display areas with enlarged pores in the skin by deepening the color of the pores. The pores in the facial skin can be seen more clearly and intuitively through the form of data. Pore clogging refers to the pores on the surface of the skin being blocked, which prevents sebum from being discharged normally, accompanied by the accumulation of stratum corneum and dirt. This phenomenon usually manifests itself in the form of blackheads, whiteheads or acne, and in severe cases may lead to skin problems such as acne and folliculitis.

Superficial Pigment

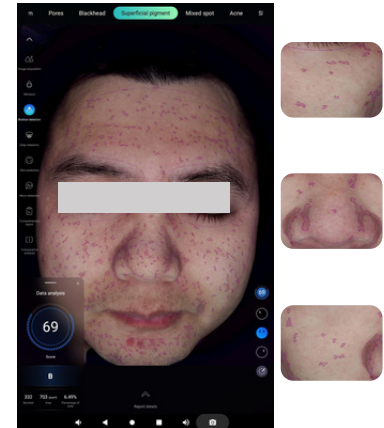


Superficial pigment refers to pigmentation that has formed on the superficial layer of the skin, including acne scars, spots, inflammatory pigmentation, etc. The coverage of pigmentation may exist in both deep and shallow layers. You can compare the image with the deep pigment image. If the shallow layer shows pigmentation but the deep layer shows no pigmentation, it means that the pigment is only deposited in the superficial layer of the skin. The algorithm marks the pigmented area with a purple polygon curve, and the shallow pigment can be seen more clearly and intuitively through the form of data.

Blackhead

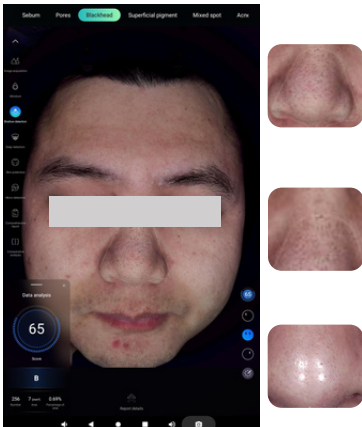


Under negative polarized light source, you can see blackheads formed by pores clogged by oil in the T-zone. The algorithm uses RBX technology to highlight the blackheads in the T-zone by deepening their color. The blackheads of the nose can be seen more clearly and intuitively through the data. Blackheads are formed by excess oil accumulation in the nose area of the skin and air oxidation. Areas with large pores are more likely to accumulate and store oil and dust in the air, so it is necessary to clean and moisturize in time to reduce the formation of large pores.

Mixed Spot

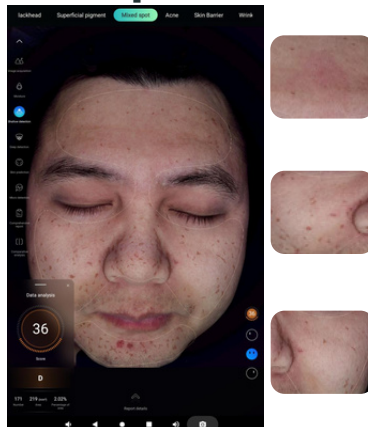


Under negative polarized light, we can see the distribution of mixed spots on the skin surface. The algorithm identifies the facial complex spot area and marks it with a brown block. The mixed spot map shows skin pigmentation such as melasma, age spots, and freckles. Melasma is a darker patch on the skin that can appear brown, black, or dark brown. Melasma may expand over time, especially if daily sun protection and skin care are not taken care of. Some melasma may be slightly raised and feel slightly convex to the touch. The algorithm marks the mixed spot area with brown color blocks, and the mixed spot situation can be seen more clearly and intuitively through the form of data.

Acne

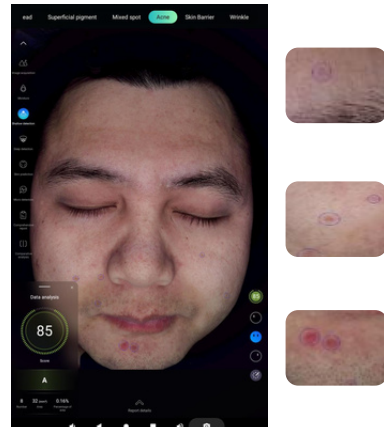


Look at the distribution of skin acne and superficial redness under negative polarized light. When pores are clogged with oil and dust, it is easy to fester or form inflammation, which will then turn into acne and acne. The algorithm identifies the distribution area of facial acne and marks it with blue circles. The more the number and the more obvious the redness of the skin, the more serious the skin acne problem is, and the skin needs to be oil-controlled to unclog the pores and eliminate inflammation.

Brown Pigment

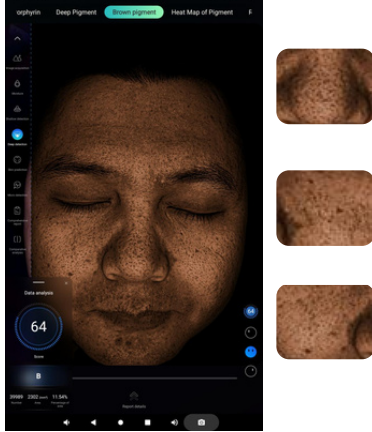


The depth of the overall brown color of the skin is mainly related to the skin color. People with darker skin or more hemoglobin have darker overall pigmentation. The areas with heavier pigmentation in the image are mostly those with higher pigment concentration density. Through the form of data, the brown pigment can be seen more clearly and intuitively.

Heat Map of Pigment

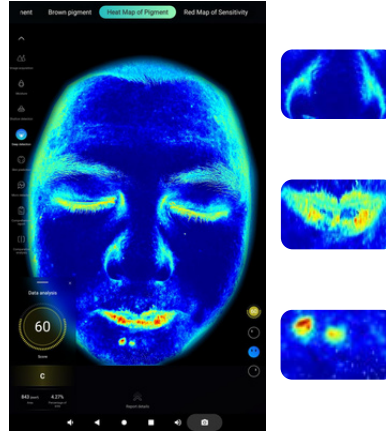


The algorithm identifies the distribution of pigments on the face and presents it in the form of a heat map. Different colors are used to represent the distribution of spots, moles, and scars visible to the naked eye under negative polarized light. Red indicates severe skin pigmentation, yellow for medium, green for lighter skin, and blue for normal skin. The situation of pigmentation can be seen more clearly and intuitively through the form of data. Pigment production mechanism: the body's own regulation, physical or chemical factors stimulate melanocytes, so that their number increases and activity increases. The melanin produced cannot be completely removed with the stratum corneum and blood circulation, and finally deposits in the local skin.

Deep Pigment

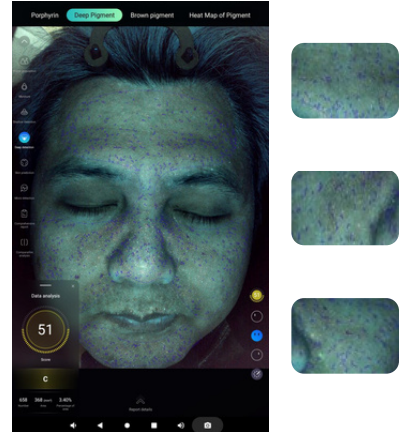


The bottom image in the picture is wood's light, and the blue-purple area is the algorithm's recognition of the facial comprehensive freckle area and the marking of the polygon curve. The situation of the deep pigment can be seen more clearly and intuitively in the form of data. The dark (black, brown) block or dot skin that appears on the face is a display of skin pigmentation (such as: melasma, freckles, malar spots, inflammatory pigmentation, acne marks, hemoglobin aggregation, etc.). The pigmentation in the deep layer of the skin can be compared with the sensitivity to determine whether it is an inflammatory hemoglobin accumulation or a spot problem.

Skin Barrier

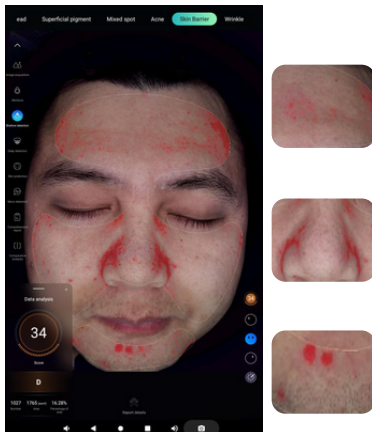


We can check the skin barrier health under negative polarized light source. The barrier image shows the skin redness problem and the distribution of red blood streaks. The formation of red blood streaks is mainly due to the damage of keratin, the weakness of the epidermis, and the long-term damage of the capillary position, which leads to vascular dilation and congestion. The red area indicates that the skin barrier is damaged, which can be used as a reference for judging the skin sensitivity and inflammation area. The damage to the barrier can be seen more clearly and intuitively through the form of data.

Wrinkle

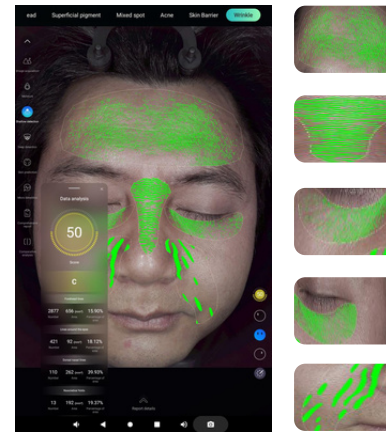


The texture of the skin surface can be viewed under a positive polarized light source. The wrinkle image shows the roughness of the skin texture, such as large pores, dry lines, fine lines, and wrinkles. It can be used as a reference for judging the fineness of the skin and the loss of collagen. The algorithm identifies the lines of the facial skin and marks the distribution of the five parts of the skin wrinkles (head-up lines, nose-back lines, peri-eye lines, crow's tail lines, decree lines) with a green short line. The more intermittent lines, the rougher the skin. You can see the wrinkles more clearly and intuitively through the form of data.

Porphyrin

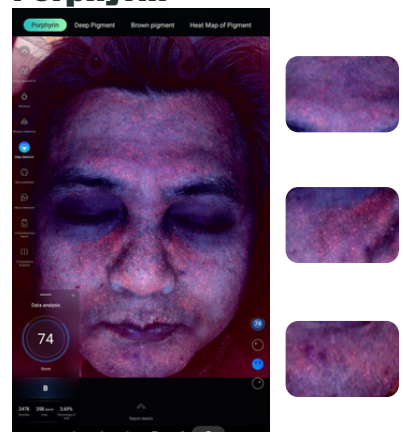


The brick-red fluorescent spots in the picture are propionibacterium acnes and malassezia. These two bacteria will aggravate the occurrence of skin acne, so they can be used as a basis for judging skin acne. Through the form of data, the situation of porin can be seen more clearly and intuitively. The living environment of propionibacterium acnes and Malassezia must have oil, so they can be used as a basis for judging the accumulation of oil in skin pores.

Red Map Sensitivity

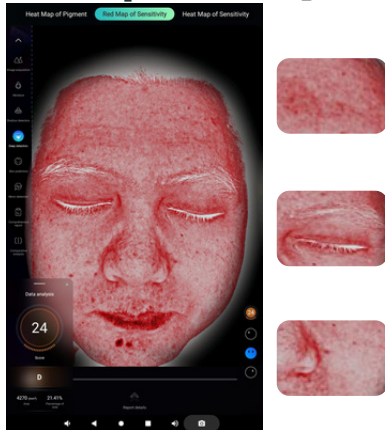


Under negative polarized light source, we can check the redness of the superficial layer of the skin and the distribution of red blood vessels. Clear distribution of red blood vessels in polarized light indicates that the skin has thin cuticle and is sensitive, requiring proper protection and care. The depth of the hemoglobin base color is related to the overall skin color. People with less hemoglobin will have a lighter color. Areas with more concentrated red represent areas where skin hemoglobin accumulation is more concentrated, which can be used as a reference for judging skin sensitivity and inflammation areas. Through the form of data, you can see the degree of skin sensitivity more clearly and intuitively.

Heat Map of Sensitivity

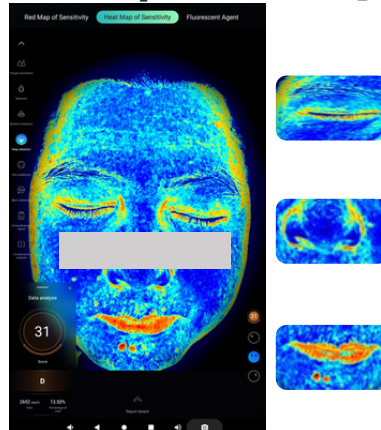


The "heat map of sensitivity" represents skin sensitivity. When the skin shows significant redness and thinning of the stratum corneum, it becomes more susceptible to external stimuli and damage, leading to issues such as dryness, sensitivity, and redness. The heat map of sensitivity is based on the distribution of subcutaneous capillaries, with areas of greater sensitivity having more capillaries. Visible redness and acne on negative polarized light images indicate areas of severe sensitivity. The algorithm uses different colors to indicate varying degrees of sensitivity and their distribution on the skin. Areas with severe sensitivity are shown in deep red, including the lips; medium sensitivity is represented in yellow, mild sensitivity in green, and normal skin appears in blue. The sensitivity is more clearly and intuitively reflected in the form of data.

Fluorescent Agent

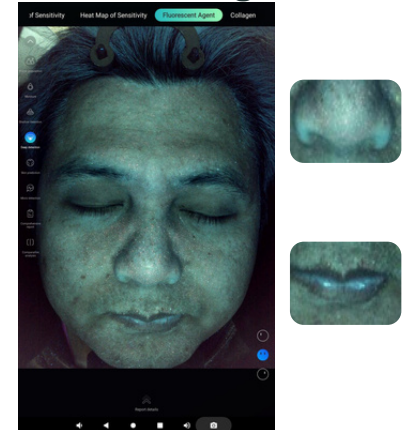





Fluorescent agent and pigments might both appear in facial imaging. To assess the fluorescent agent, focus specifically on the fluorescence response. The difference between fluorescent agent and porphyrins is as follows: porphyrins exhibit brick-red fluorescent agent, while fluorescent agent display intense blue light and usually appear as large, sheet-like areas. The difference between fluorescent agent and facial dust is as follows: facial dust appears as white, bright, floating, and short, wispy lines on the surface, while fluorescent agent typically display bright colors and are often more diffuse or spread over larger areas.

Collagen

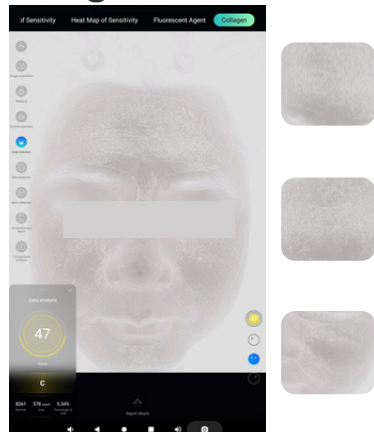

Under polarized light, we can assess the loss of collagen on the skin's surface. Mixed light images reveal skin texture issues such as enlarged pores, dry lines, fine lines, and wrinkles. They serve as a reference for evaluating skin smoothness and collagen loss. In mixed light images, a higher number of discontinuous lines indicates rougher skin texture and more severe collagen loss. The sensitive process is more clearly and intuitively reflected in the form of data.

Moisture test

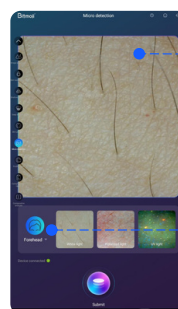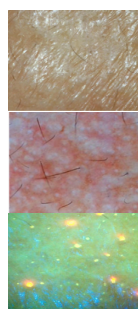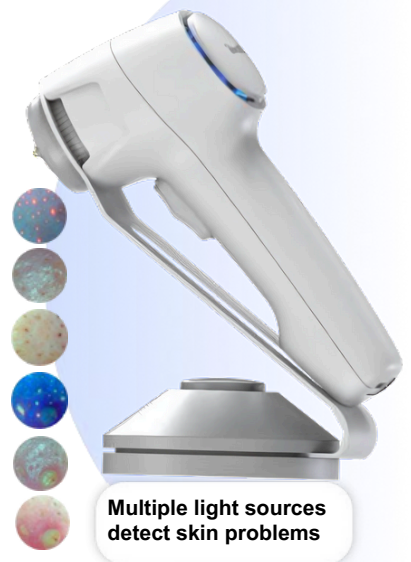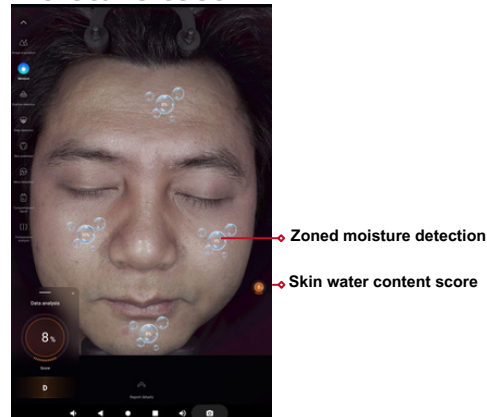

Local microscopic display

Location selection, partition detection.



Detection function-Moisture test report 2

Details of the report

Problem analysis

Your skin is in a state of lack of water, sebum secretion is less, the skin is not shiny, dry and desquamate, prone to small wrinkles and stains, poor tolerance to external stimuli

Natural moisturizing factors and sebum secrete less in the keratinocytes, the skin is thinner, dry and lackluster, and easily dandruff. After washing the face, it often feels tight, sensitive to external stimuli such as sun exposure and wind blowing, lacks elasticity, and is prone to aging, wrinkles and pigmentation. Dry skin is not only related to congenital factors, but also related to external factors such as frequent exposure to the sun and wind, excessive use of alkaline detergents, abuse of cosmetics, etc. Improper care is easy to turn sensitive skin.

Daily Care Advice

1. Diet: daily drinking water is recommended to maintain 2000ml-2500ml, eat more vegetables, fruits and vegetables and gluten-rich foods, to supplement water.
2. Daily care: according to their own needs to choose appropriate skin care products. Emulsion water, suitable for oily skin or summer use; The cream is oily and thick, suitable for dry skin or winter use.
3. Healthy work and rest Ensure adequate sleep and regular work and rest, avoid the body endocrine disorders and thus affect the health of the skin and moisture.

Program Recommendation

Dermapen

A microneedling device that boosts collagen to improve scars, wrinkles, and skin texture.



Mesotherapy for pigmentation:

Targets dark spots with skin-brightening ingredients to even out skin tone.



HydraFacial:

Deeply cleanses, exfoliates, and hydrates for instantly glowing skin.



Toner:

Balances skin pH and refreshes, prepping it for better absorption of skincare.



Botox and fillers for wrinkles:

Botox relaxes lines; fillers add volume to smooth deeper wrinkles.



Laser types:

Ablative, non-ablative, and pigment/vascular lasers for various skin concerns.




scan the QR code to view the report in your mobile



Key Features of the EPILYUM AI REVEAL :

- Uses 8 advanced spectral imaging technologies
- Powered by intelligent AI for one-click analysis
- Detects 17+ skin conditions like wrinkles, acne, and pigmentation
- Delivers clear, data-driven reports with 3D imaging and heat maps
- Supports cloud storage

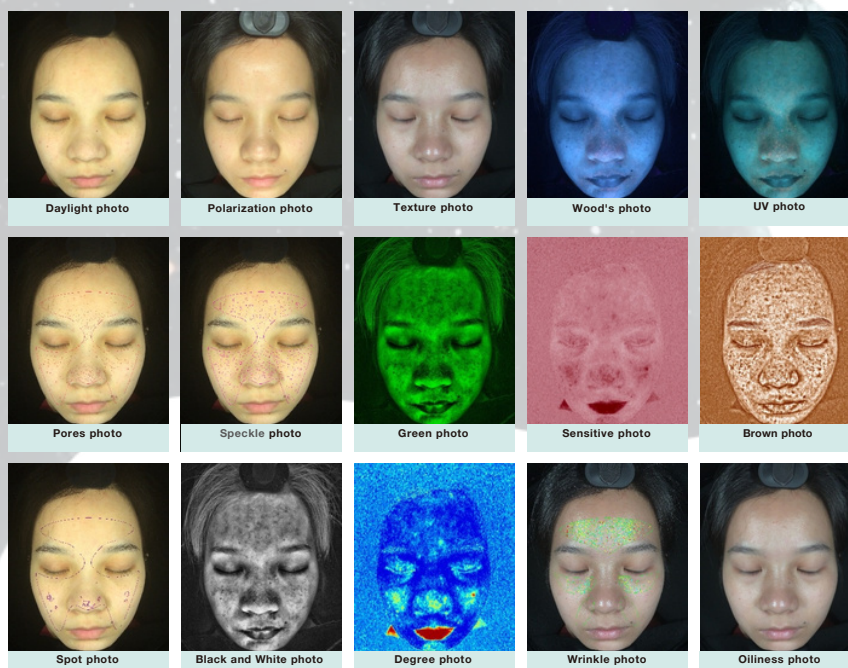


Skin prediction




Deeply predict the future of the skin and awaken patient's desire for young skin.

According to the current skin condition of the customer, through the training of AI large model, simulate the skin condition after customer care and the aging situation of different age groups.

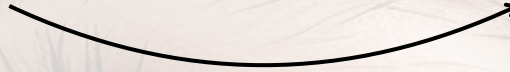

Cloud storage file management

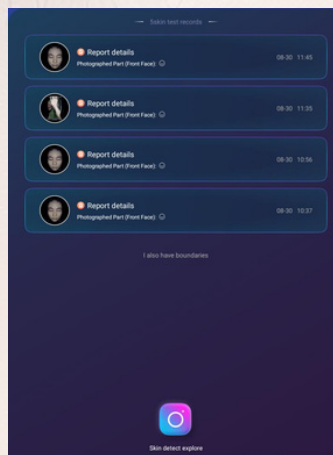

- File Management
- Profile
- One click search
- Cloud storage
- Number of skin tests

Instrument Specifications

EPILYUM - AI SKIN ANALYZER	Input Power: 24V--2A	Power Consumption: 40W
Lighting Modes: RGB (Standard Daylight), Cross-Polarized Light, Parallel Polarized Light, UV (Ultraviolet Light), Woods Light		
Lighting Technology: Solid State LED	Operating Environment: 5~40°C	Weight: 15kg

