



ThermiQ Microneedling RF

Microneedling + Multipolar RF
diversified treatment options , excellent results



EPILYUM ThermiQ Microneedling RF

ThermiQ Microneedling RF is an advanced skin rejuvenation system that combines precise microneedling with radiofrequency energy. It works by creating microchannels in the skin while delivering controlled heat deep into the dermis, stimulating collagen and elastin production.

✓ Advanced Skin Rejuvenation:

Combining microneedling with radiofrequency energy, **ThermiQ** stimulates collagen and elastin production to reduce wrinkles, tighten skin, improve elasticity, and minimize scars, pores, and stretch marks—all with minimal downtime.

✓ Precision & Safety:

Fine microneedles deliver controlled RF energy deep into the dermis while protecting the epidermis, ensuring effective treatment with little discomfort, rapid healing, and no risk of pigmentation. Adjustable depth, energy, and speed settings allow tailored solutions for different skin types and conditions.

✓ Comprehensive Skin Solutions:

From wrinkle reduction and skin lifting to acne scar repair and stretch mark improvement, **ThermiQ** offers a full range of non-surgical rejuvenation treatments, making it a versatile solution for both facial and body care.

✓ Microneedle handle with vacuum suction (negative pressure handle)

The handpiece incorporates advanced vacuum-assisted (negative pressure) adsorption technology to stabilize the skin surface during treatment. This mechanism gently lifts the skin before needle insertion, ensuring uniform needle penetration depth and precise RF energy delivery across the treated area. The system minimizes patient discomfort, enhances treatment consistency, prevents drug leakage, and allows for automatic needle retraction after each pulse. This innovative design ensures safer, more effective, and scientifically controlled treatment outcomes.



Multiple Needle Options - Range: 0.1 - 4.0mm

10 Pins: Designed for small, delicate areas (around eyes, nose, or fine wrinkles).

25 Pins: Ideal for medium areas such as cheeks, forehead, and general facial rejuvenation.

64 Pins: Best for larger treatment zones like the neck, body, or areas requiring deeper collagen stimulation.



10 pins



25 pins




64 pins

✓ Microneedle Handle without Vacuum Suction

Efficient & Minimal Trauma:

The non-vacuum handle uses ultra-fine microneedles (0.23mm) to create precise micro-injuries with minimal discomfort and rapid healing.

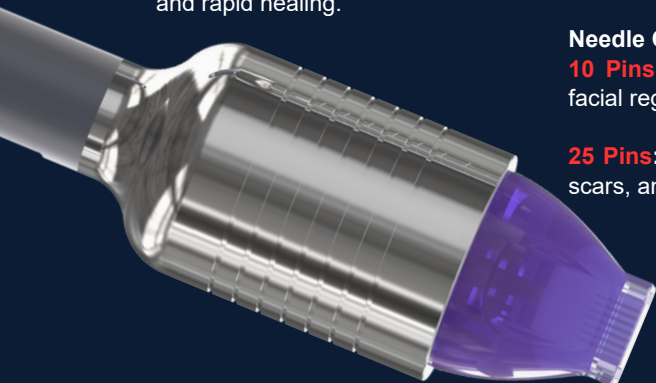

Needle Options - Range: 0.1 - 8.0mm

10 Pins: Suitable for smaller or targeted areas such as fine lines and delicate facial regions.

25 Pins: Designed for medium coverage on the face, helping with wrinkles, acne scars, and overall skin rejuvenation.



10 pins



25 pins

Deep Precision Technology for Exceptional Skin Tightening

8MM

In our new device, the needle penetration depth reaches up to **8 mm**, allowing precise delivery of energy into deeper dermal layers. This advanced capability significantly enhances collagen remodeling and tissue tightening, resulting in superior clinical outcomes and remarkable skin-firming effects.




smaller wound:

With a needle diameter of only **0.23 mm**, the system minimizes tissue damage during penetration, resulting in significantly smaller **micro-wounds**. This advanced design promotes faster healing, reduces post-treatment downtime, and ensures a more comfortable recovery experience for patients.



insulation coating :

Each needle column is precisely engineered with a **specialized insulation coating** that serves a dual purpose: it concentrates the radiofrequency (RF) energy exclusively at the needle tip while shielding the surrounding skin surface from unnecessary thermal exposure. This advanced insulation design allows the RF energy to penetrate deeply and accurately into the targeted tissue layers, enhancing collagen remodeling and skin tightening effects. At the same time, it prevents surface burns or epidermal damage, ensuring a safer, more comfortable, and highly effective treatment experience for patients.

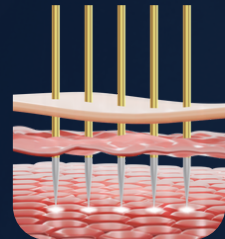

Some exemplary treatment outcomes achieved with the **ThermiQ Microneedling RF, along with their optimized parameters:**



PULSE WIDTH 0.4ms
NEEDLE LENGTH 1.8mm
SPEED 4ms
ENERGY 30mj
3 SESSIONS , 2 TIMES



PULSE WIDTH 0.6ms
NEEDLE LENGTH 2.4mm
SPEED 0.6ms
ENERGY 40mj
3 SESSIONS , 2 TIMES




PULSE WIDTH 0.3ms
NEEDLE LENGTH 2mm
SPEED 5ms
ENERGY 40mj
5 SESSIONS , 2 TIMES




PULSE WIDTH 0.4ms
NEEDLE LENGTH 1.8mm
SPEED 5ms
ENERGY 40mj
3 SESSIONS , 2 TIMES


Imported
accessories

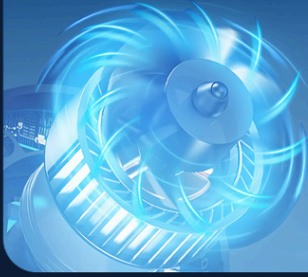


Operable in multiple
parts

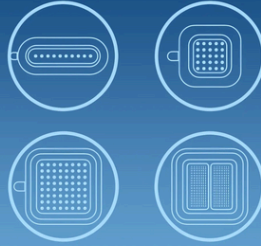

Safe and
reliable

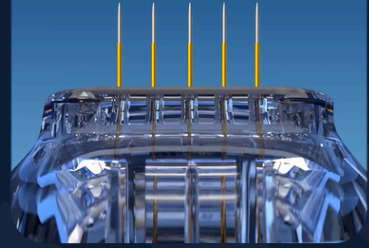

Negative pressure adsorption technology

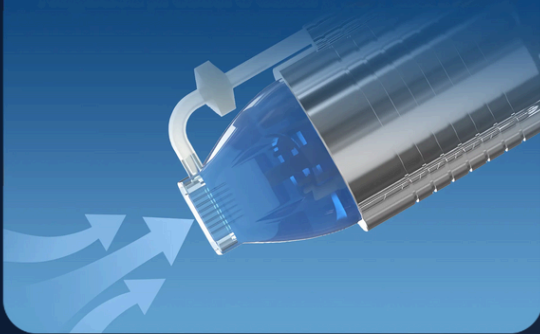

Fast
treatment



TECHNICAL FEATURES

Screen	12.1 inch LCD
Size	61*57*120cm
RF Frequency	2MHZ
Speed	1-10 level optional
Pressure Level	1-10 level optional
Needle Depth	(0.1 - 8.0) mm adjustable
Needle Coating	Insulated
Tips	M-0,M-10,M-25,M-64
Working mode	Automatic /Manual
Handle	Negative Pressure Handle / Quick Needle Handle

Warnings and risks : RF Microneedling devices are intended solely for use by qualified physicians or practitioners trained in the use of radiofrequency microneedling technology. Incorrect treatment settings or misuse of the device may present a risk of serious injury to both the patient and the operator. Potential risks associated with any RF Microneedling procedure may include changes in pigmentation, infection, redness (erythema), swelling (induration), or scarring. Always read and fully understand the RF Microneedling system and accessory operator manuals for the complete list of indications, contraindications, and possible risks.

